

C400

Four-channel Pulse Counting Detector Controller

User Manual



Four-channel Pulse Counting Detector Controller

CONTENTS

FIGURES	6
SAFETY INFORMATION	7
MODELS	9
SCOPE OF SUPPLY	10
OPTIONAL ITEMS	11
POWER SUPPLIES	11
SIGNAL CABLES AND CABLE ACCESSORIES	11
DATA CABLES.....	11
POWER AND CONTROL CABLES.....	11
PRE-AMPLIFIERS	11
FIBER-OPTIC LOOP.....	12
INTENDED USE AND KEY FEATURES.....	13
INTENDED USE.....	13
KEY FEATURES	13
SPECIFICATION.....	14
INSTALLATION.....	19
MOUNTING	19
GROUNDING AND POWER SUPPLY	19
CONNECTION TO SIGNAL SOURCE	20
<i>Typical setup</i>	20
<i>Cables</i>	22
HOW THE C400 WORKS - AN OVERVIEW	23
PULSE HEIGHT DISCRIMINATION AND COUNTING	23
TRIGGERING.....	23
ACQUISITION MODES	24
DETECTOR SERVICES	24
SELF-TESTING AND CALIBRATION	24
MONITOR OUTPUTS	25
COMMUNICATION TO THE HOST COMPUTER	25
EMBEDDED SOFTWARE	25
GETTING STARTED USING THE PTC DIAGNOSTIC G2 HOST PROGRAM.....	26
INSTALLATION.....	26
CONNECTING TO THE C400	27
SCREEN LAYOUT.....	31
GRAPHIC DISPLAY.....	32

Four-channel Pulse Counting Detector Controller

DATA LOGGING.....	35
DATA DISPLAY.....	35
SETUP DISPLAY: DETECTOR, DISCRIMINATOR AND SCALER CONTROL	36
SETUP 2 DISPLAY: TEST PULSER AND ENCODER INPUT CONTROLS	38
DISCRIMINATOR SWEEP.....	39
CALIBRATION	40
PROPERTIES.....	41
<i>Firmware update.....</i>	41
<i>Changing the IP address</i>	43
<i>Restoring the default IP address.....</i>	43
C400 CIRCUIT OVERVIEW	44
DISCRIMINATOR INPUTS	44
OVERALL CIRCUIT	44
HIGH VOLTAGE SUPPLIES.....	47
SETTING THE HIGH VOLTAGE SUPPLIES	47
CHANGING THE HIGH VOLTAGE SUPPLY RANGE AND POLARITY	48
TRIGGERS AND DATA BUFFERING.....	49
C400 ACQUISITION STATES.....	49
TRIGGER MODES.....	50
BUFFER MEMORY.....	51
MONITOR OUTPUTS.....	53
DISCRIMINATOR SWEEP	54
DEADTIME CORRECTION.....	56
MODELS OF DEADTIME	56
DETERMINING THE DEADTIME	57
PRE-AMPLIFIER SELECTION	59
PULSER OUTPUTS.....	60
ENCODER INPUT.....	61
CONNECTORS	63
FRONT PANEL CONNECTORS	63
<i>Signal input</i>	63
<i>Pre-amplifier power and other services.....</i>	63
<i>High voltage outputs</i>	63
REAR PANEL CONNECTORS.....	63
<i>Ethernet communications.....</i>	63
<i>Fiber-optic communications:.....</i>	64
<i>RS-232 / RS-485 communications.....</i>	64
<i>Gate output</i>	64

Four-channel Pulse Counting Detector Controller

<i>Gate input</i>	64
<i>Encoder input</i>	65
<i>Monitor outputs</i>	65
<i>Power input</i>	65
<i>Ground lug</i>	65
CONTROLS AND INDICATORS	66
FRONT PANEL CONTROLS.....	66
REAR PANEL CONTROLS	66
<i>Reset button</i>	66
<i>Address switch</i>	66
<i>Mode switch</i>	66
FRONT PANEL INDICATORS	67
<i>HV on</i>	67
REAR PANEL INDICATORS	67
24V	67
2.5V	67
<i>Initiated</i>	67
<i>Active</i>	67
<i>Com</i>	67
<i>Ethernet</i>	67
<i>Optical</i>	67
<i>Serial</i>	67
<i>Power up sequence</i>	68
<i>RJ-45 indicators</i>	68
COMMUNICATIONS INTERFACES	69
RS-232 CONNECTION	69
RS-485 CONNECTION EXAMPLE	70
ASCII COMMUNICATION	71
IEEE 488.2 MANDATORY COMMANDS	72
IEEE 488.2 OPTIONAL COMMANDS.....	72
C400 SPECIFIC COMMANDS	73
DIGITAL STATUS BITS	75
USING PUTTY TO TEST THE SERIAL ASCII CONNECTION	75
CREATING A HOST SOFTWARE APPLICATION	77
FAULT-FINDING	78
MAINTENANCE	82
RETURNS PROCEDURE	83
SUPPORT	84

Four-channel Pulse Counting Detector Controller

DISPOSAL 85

DECLARATION OF CONFORMITY 86

REVISION HISTORY..... 87

Four-channel Pulse Counting Detector Controller

FIGURES

Figure 1	Serial adaptor cable	11
Figure 2	CP10 pre-amplifier	12
Figure 3	C400 front and rear panels. Dimensions mm.	17
Figure 4	C400 case plan and side views. Dimensions mm.	18
Figure 5	Schematic example C400 installation	21
Figure 6	Comparator operation	23
Figure 7	PTC DiagnosticG2 installation.....	27
Figure 8	Direct Ethernet connection.....	28
Figure 8a	Configuring the IP address on the host computer.....	28
Figure 9	Ping test of the Ethernet connection.....	29
Figure 10	Discover devices	29
Figure 11	Ready to connect to the C400.....	30
Figure 12	Opening the C400 window	30
Figure 13	Strip chart data display showing cursor. The scope display looks similar.	32
Figure 14	Histogram data display showing cursor.....	33
Figure 15	Increasing the low pass filtering of the displayed data.	34
Figure 16	Zero offset subtraction inactive and active.	34
Figure 17	Data display.....	35
Figure 18	Primary setup parameters.....	36
Figure 19	Secondary setup parameters.....	38
Figure 20	Discriminator sweep controls.....	39
Figure 21	Calibration control and display.....	40
Figure 22	Firmware update warning.	41
Figure 23	Selecting the firmware file.	42
Figure 24	Firmware update status window.	42
Figure 25	IP address configuration dialog.....	43
Figure 26	The basic discriminator circuit.....	44
Figure 27	C400 block schematic.....	45
Figure 28	High voltage jumper settings	48
Figure 29	C400 acquisition state diagram.	49
Figure 30	Unbuffered logged data with a 10 Hz comms rate, 100000 μ sec integration time.....	51
Figure 31	Unbuffered logged data with a 10 Hz comms rate, 10000 μ sec integration time.....	51
Figure 32	Buffered logged data with a 10 Hz comms rate, 10 μ sec integration time.....	52
Figure 33	Example of monitor output pulse.....	53
Figure 34	Na-22 spectrum acquire using discriminator sweep.....	54
Figure 35	Non-paralysable deadtime model.	56
Figure 36	Deadtime measurement using a controllable intensity source.	57
Figure 37	Pulse shapes from a LaCl ₃ scintillation detector as seen by the C400 with CP10-A and CP10-B pre-amplifiers.....	59
Figure 38	Test pulser LED connection.	60
Figure 39	Differential quadrature signals from a shaft encoder, clockwise movement.	61
Figure 40	Encoder signals when no differential signals are available.....	62
Figure 41	Example bias circuit for encoder inputs when no differential signals are available.	62
Figure 42	RS-232 cable C400 to PC.....	69
Figure 43	RS-485 cable C400 to TC100.	70
Figure 44	Setup of PuTTY for ASCII serial communication with the C400.	76
Figure 45	Example terminal session.....	76

Four-channel Pulse Counting Detector Controller

SAFETY INFORMATION

This unit is designed for compliance with harmonized electrical safety standard EN61010-1:2000. It must be used in accordance with its specifications and operating instructions. Operators of the unit are expected to be qualified personnel who are aware of electrical safety issues. The customer's Responsible Body, as defined in the standard, must ensure that operators are provided with the appropriate equipment and training.

The unit is designed to make measurements in **Measurement Category I** as defined in the standard.



CAUTION. According to installed options the C400 can generate high voltages as follows, present on the central conductor of the SHV (Safe High Voltage) connectors:

+ or – 2000 V DC at 0.5 mA maximum.

or

+ or – 1000 V DC at 1.0 mA maximum

or

+ or – 500 V DC at 2.0 mA maximum

or

+ or – 200 V DC at 5.0 mA maximum

The hazardous live voltages on the SHV central conductor are not accessible under the definitions of EN61010 but may nevertheless give a noticeable shock if misuse were to lead you to come into contact with them. The user must therefore exercise appropriate caution when using the device and when connecting cables. Power should be turned off before making any connections.

The unit must not be operated unless correctly assembled in its case. Protection from high voltages generated by the device will be impaired if the unit is operated without its case. Only Service Personnel, as defined in EN61010-1, should attempt to work on the disassembled unit, and then only under specific instruction from Pyramid Technical Consultants, Inc..

The unit is designed to operate from +24VDC power, with a maximum current requirement of 1000 mA. A suitably rated power supply module is available as an option. Users who make their own power provision should ensure that the supply cannot source more than 3000 mA.

A safety ground must be securely connected to the ground lug on the case.

Some of the following symbols may be displayed on the unit, and have the indicated meanings.

Four-channel Pulse Counting Detector Controller



Direct current



Earth (ground) terminal



Protective conductor terminal



Frame or chassis terminal



Equipotentiality



Supply ON



Supply OFF



CAUTION – RISK OF ELECTRIC SHOCK



CAUTION – RISK OF DANGER – REFER TO MANUAL

Four-channel Pulse Counting Detector Controller

MODELS

C400	Four-channel pulse-counting detector controller
-nXP20/10/5/2	Add positive 0 to 2000V / 1000 V / 500 V / 200 V HV outputs to n channels
-nXN20/10/5/2	Add negative 0 to 2000V / 1000V / 500 V / 200 V HV outputs to n channels

Example:

C400-4XP10

C400 with four 1000V positive auxiliary high voltage outputs.

Four-channel Pulse Counting Detector Controller

SCOPE OF SUPPLY

C400 model as specified in your order.

PSU24-40-1 24 VDC power supply.

USB memory stick containing:

- C400 User manual

- PTC DiagnosticG2 software installation files

- Terminal emulation program for testing ASCII communication

- Test data

Optional items as specified in your order.

Four-channel Pulse Counting Detector Controller

OPTIONAL ITEMS

POWER SUPPLIES

PSU24-45-1 +24 VDC 45W PSU (universal voltage input, plug receptacle for standard IEC C14 three-pin socket) with output lead terminated in 2.1mm threaded jack.

PD-8 Eight output +24 VDC power supply unit, 19" rack mounting

SIGNAL CABLES AND CABLE ACCESSORIES

CAB-LCO-10-LCO Cable, coaxial signal, Lemo 00 plug to Lemo 00 plug, 10 m.

CAB-SHV-10-SHV Cable, coaxial HV, SHV to SHV, 10 m.

DATA CABLES

AB450K-R RS-232 6 pin DIN male to 9 pin D sub female adaptor.

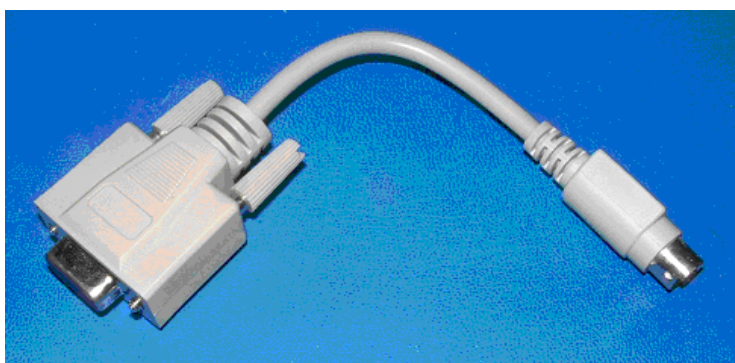


Figure 1. Serial adaptor cable

Fiber-optic cable pair 1 mm plastic fiber ST terminated with color-coded sleeves, 5 m. Other lengths available to order.

Fiber-optic cable pair 200 μ m silica fiber ST terminated with color-coded sleeves, 5 m. Other lengths available to order.

POWER AND CONTROL CABLES

CAB-L4-10-D9 Cable, screened, four pin Lemo 0B plug to D9 male, 10 m, for pre-amp power.

PRE-AMPLIFIERS

CP10-A Wide-band DC-coupled preamplifier suitable for NaI, LaCl₃, LaBr₃, and YAP scintillators coupled to photomultiplier tubes.

CP10-B Very wide-band DC-coupled preamplifier suitable for plastic scintillators, continuous dynode electron multipliers and other pulse counting detectors with very short output pulses.

Four-channel Pulse Counting Detector Controller

FIBER-OPTIC LOOP

A200 USB to fiber-optic adaptor.

A300 fiber-optic loop controller / Ethernet adaptor.

A500 intelligent real-time cell controller for five loops with Ethernet interface.

Four-channel Pulse Counting Detector Controller

INTENDED USE AND KEY FEATURES

INTENDED USE

The C400 is intended to count fast pulses generated by photomultipliers, avalanche photodiodes and similar devices on four parallel channels. It provides high voltage for the sensors and low voltage power for pre-amplifiers. Multiple C400 units can measure synchronously using a distributed trigger, and encoder pulses can be measured to coordinate pulse counting with the position of a motion system.

The C400 has design features which make it tolerant of electrically noisy environments, but the place of use is otherwise assumed to be clean and sheltered, for example a laboratory or light industrial environment. The unit may be used stand-alone, or networked with other devices and integrated into a larger system. Users are assumed to be experienced in the general use of precision electronic circuits for sensitive measurements, and to be aware of the dangers that can arise in high-voltage circuits.

KEY FEATURES

Four fully-parallel high speed window comparators and 100 MHz scalars.

Four independent high voltage outputs for detector biasing

Four power supply outlets for remote pre-amplifiers.

Four analog voltage outputs for remote HV control or other uses.

Four LED pulser outputs.

Input for quadrature position encoder.

Discriminated monitor pulse outputs.

Gate inputs and output for triggering on external events, and trigger distribution.

Analog monitor output which can track any arithmetic sum,

Ethernet interface with TCP/IP and UDP messaging.

Alternative RS-232, RS-485 and fast fiber-optic serial interfaces built-in.
Selectable baud rates.

Can be operated in a fiber-optic serial communication loop with up to fourteen other devices.

Four-channel Pulse Counting Detector Controller

SPECIFICATION

<i>Signal inputs</i>		
	Number	Four, independent parallel
	Input impedance	50 ohm
	Pulse height range	0 to 5 volt, either polarity
	Minimum pulse width	Less than 10 nsec
	Pulse pair resolution	Better than 10 nsec
	Discriminator analog bandwidth	DC - 150 MHz
	Discriminator range	Window discriminator, upper and lower levels adjustable 0 to 5 V, either polarity 16 bit resolution over bipolar range
<i>Scalers</i>		
	Number	Four
	Depth	32 bit (transparent background readout)
	Max counting rate	100 MHz
	Integration period	Configurable, 100 µsec to 1000 sec.
	Deadtime correction	Real-time non-paralyzable deadtime correction can be applied to each individual integration period.
	Internal buffering	On-board data buffer for up to 65536 integrations
<i>Test pulsers</i>		
	Number	Four
	Voltage	0 - 5 V
	Current capability	30 mA
	Pulse duration	1 nsec to 1 µsec
	Frequency	1 Hz to 1 MHz
<i>Preamplifier outputs</i>		
	Number	Four
	Power output	+12 VDC 100 mA, -12 VDC 100 mA, fused

Four-channel Pulse Counting Detector Controller

	Analog output	-5 to +5 VDC 16 bit resolution for remote high voltage control or other purpose
	Test pulser	TTL output with controllable frequency, pulse width and pulse count for test pulser LEDs
<i>High voltage outputs</i>		
	Number	Up to four, independently configurable for range and polarity at time of purchase.
	Output power	1 watt
	Voltage options	20 to 200 V. Line <0.01%, Load <0.05%, Ripple <0.01% 50 to 500 V. Line <0.01%, Load <0.01%, Ripple <0.01% 100 to 1250 V. Line <0.001%, Load <0.005%, Ripple <0.001% 200 to 2000 V Line <0.01%, Load <0.01%, Ripple <0.001%
<i>Monitor outputs</i>		
	Number	Four
	Levels	TTL levels into minimum impedance of 50 ohm
	Pulse width	10 nsec minimum
	Latency	Less than 80 nsec after discriminator conditions are met
<i>Triggers and gating</i>		
	Trigger modes	Internal (autorun), external start, external start/stop, external start/hold, external windowed
	Gate input	TTL levels into 10 kohm
	Gate output (passthrough)	TTL levels into minimum impedance of 50 ohm
	Latency	Less than 50 nsec between gate input and gate output
<i>Encoder input</i>		
	Signals	Differential position pulses (A+/A-, B+/B-) and index (Z+/Z-)
	Levels	0.2 V to 20 V differential
	Isolation	8 kV
	Maximum count rate	1 MHz
	Counter depth	32 bit (transparent background readout)
<i>Communication interfaces</i>		
	Ethernet	Ethernet 10/100 Base T, TCP/IP and UDP

Four-channel Pulse Counting Detector Controller

	RS-232	115.2 kbps, serial ASCII protocol
	RS-485	115.2 kbps, serial ASCII protocol
	Fiber optic	10 Mbps, binary serial protocol. For integrating into Pyramid systems architectures.
<i>Controls and indicators</i>		
	Controls	Rotary switches for communications mode and fiber optic loop address, rear panel Push button processor reset, rear panel
	Indicators	Four HV enabled yellow LEDs, front panel Eight status green LEDs, rear panel
<i>Power input</i>		
	Voltage	+24 VDC (+/-2 V)
	Current	Up to 1000 mA
<i>Case</i>		
	Format and materials	1U 19" 250 mm deep rack mounting steel chassis with Al alloy front panel and polycarbonate decals. Fan-cooled.
	Protection rating	The case is designed to rating IP43 (protected against solid objects greater than 1mm in size, protected against spraying water).
	Weight	3.5 kg (7.7 lb)
	Dimensions	(see figures 1 and 2).
<i>Environment</i>		
	Operating	0 to 35 C (15 to 25 C recommended) < 70% humidity, non-condensing vibration < 0.1g all axes (1 to 1000Hz)
	Shipping and storage	-10 to 50C < 80% humidity, non-condensing vibration < 2 g all axes, 1 to 1000 Hz

Four-channel Pulse Counting Detector Controller

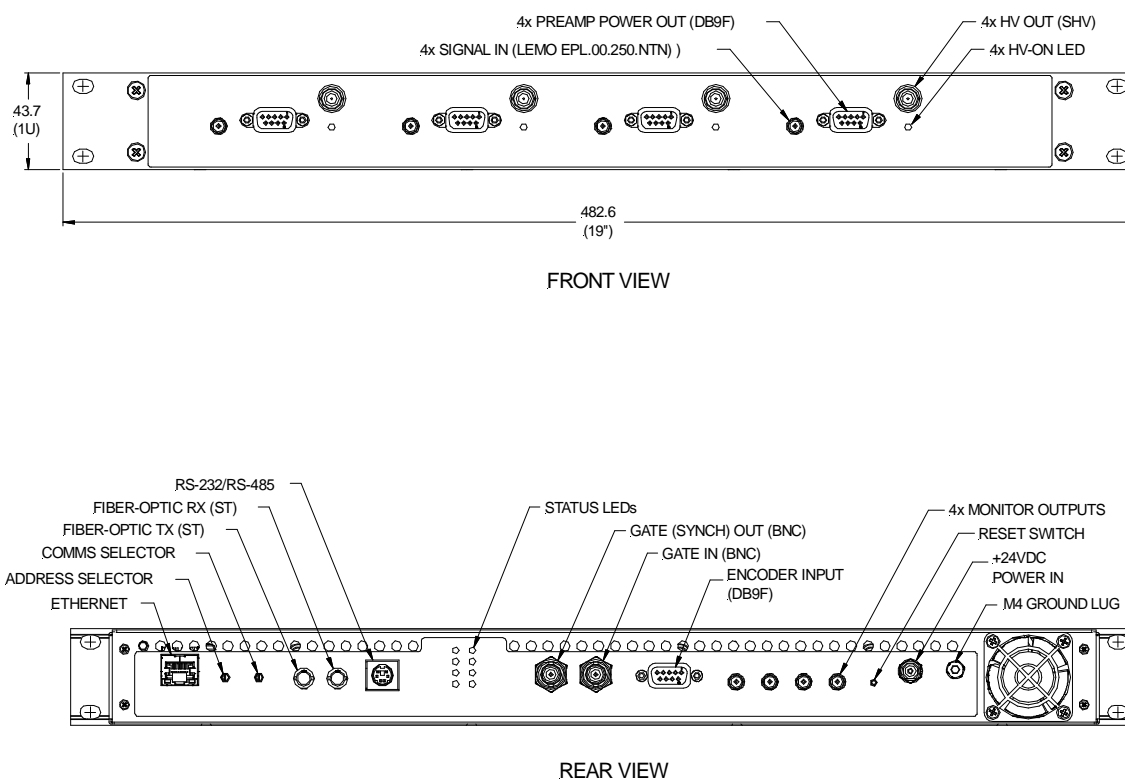


Figure 3. C400 front and rear panels. Dimensions mm.

Four-channel Pulse Counting Detector Controller

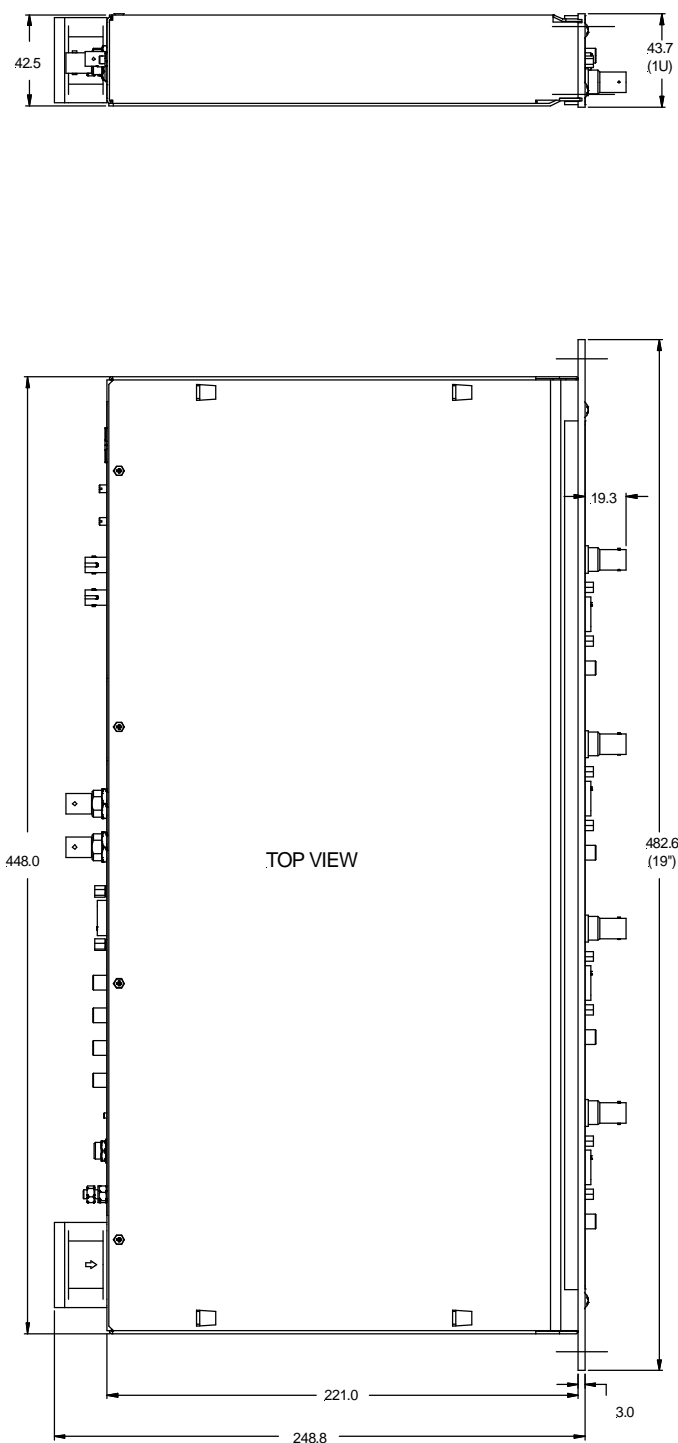


Figure 4. C400 case plan and side views. Dimensions mm.

Four-channel Pulse Counting Detector Controller

INSTALLATION

MOUNTING

The C400 is intended for 19" rack mounting, but may be mounted in any orientation, or may be simply placed on a level surface. A fixed mounting to a secure frame is recommended in a permanent installation for best low current performance, as this can be degraded by movement and vibration.

The mounting position should allow sufficient access to connectors and cable bend radii. 60 mm minimum clearance is recommended at front and back of the device.

Best performance will be achieved if the C400 is in a temperature-controlled environment. No forced-air cooling is required in addition to the unit's built-in fan, but free convection should be allowed around the back and sides of the case.

GROUNDING AND POWER SUPPLY

A secure connection should be made using a ring lug, from the M4 ground lug to local chassis potential. This is the return path for any high voltage discharge passing via the C400.

+24 VDC power should be provided from a suitably-rated power supply with the following minimum performance:

Output voltage	+24 +/- 0.5 VDC
Output current	1000 mA minimum, 3000 mA maximum
Ripple and noise	< 100 mV pk-pk, 1 Hz to 1 MHz
Line regulation	< 240 mV

The C400 includes an internal automatically re-setting PTC fuse rated at 1.1 A. However the external supply should in no circumstances be rated higher than the connector limit of 5 A.

Four-channel Pulse Counting Detector Controller

CONNECTION TO SIGNAL SOURCE

Typical setup

Figure 3 shows an example installation in schematic form. A scintillator / photomultiplier generates fast pulses when ionizing radiation interacts with the scintillator. The pulses are negative going and ground referenced, because we are measuring a current of electrons, and the photomultiplier is biased with negative high voltage at the photocathode. A pre-amplifier such as the Pyramid CP10 amplifies the pulses and sends them over a coaxial cable to the C400 input. Here the pulse height is compared to the lower and upper discriminator levels, and if it lies between them, a pulse is counted. High voltage for the photomultiplier is provided by the C400, as is DC power for the pre-amplifier. A test LED built into the scintillator assembly is driven by the pulser output of the C400 to provide a diagnostic test.

In the example shown, a goniometer is used to position a sample, for example in a diffraction experiment, and the C400 is able to read the position encoder so that signal counts can be correlated to sample position. A gate input signal is used to trigger measurements, and the C400 gate output signal is passed to another unit that needs to be triggered, for example another C400. The monitor output is used to feed pulses to an external counting device, such as a ratemeter.

Four-channel Pulse Counting Detector Controller

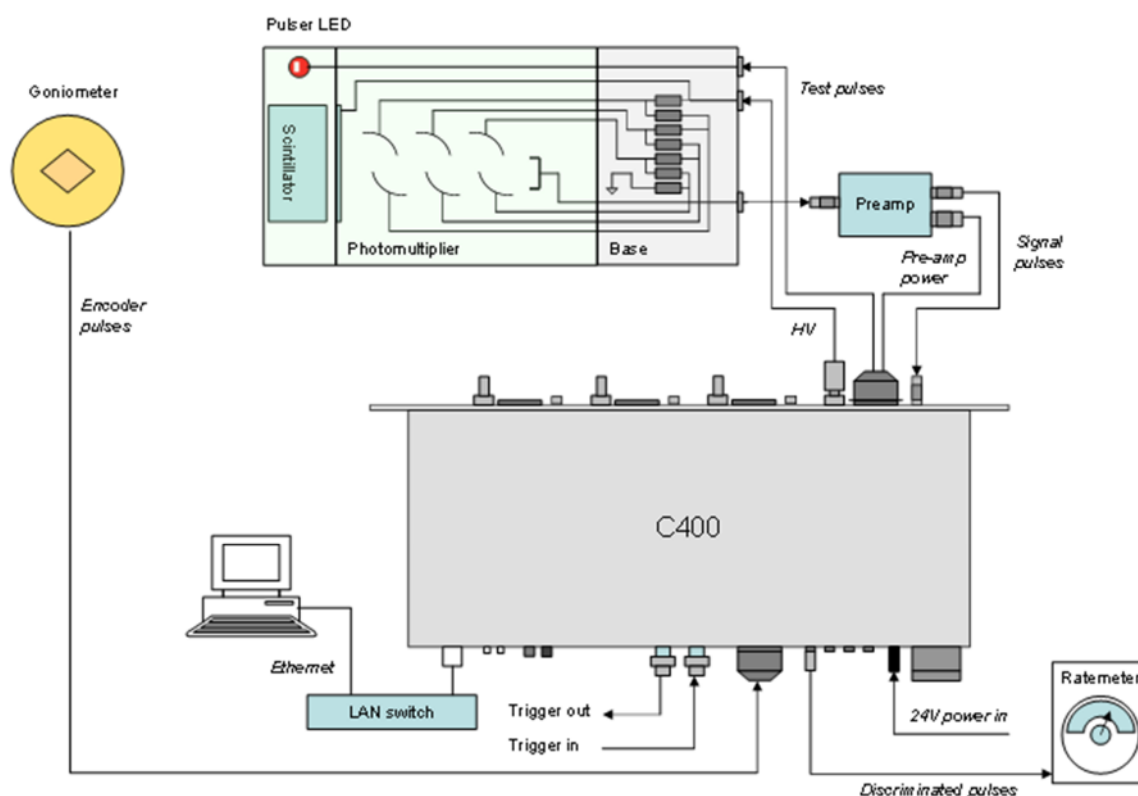


Figure 5. Schematic example C400 installation

Other signal sources can include channel electron multipliers, avalanche photodiodes, silicon photomultipliers and diamond detectors, or indeed any sensor that produces fast pulses in response to incoming radiation of particles or photons.

Many photomultipliers are biased with negative high voltage as shown in figure 3, with the result that the pulses emerging at the bottom of the dynode chain are ground-referenced, and can be connected directly to DC-coupled pre-amplifier, or even directly to the C400 inputs if they have sufficient amplitude. However, some photomultipliers are biased with positive high voltage at the anode and a grounded photocathode. The pulses are then referenced to the high voltage. In this case they must be AC-coupled using a suitably-rated capacitor before connecting to the preamplifier or the C400. Often such decoupling capacitors are integrated into the tube base bias network.

Four-channel Pulse Counting Detector Controller

Cables

A pre-amplifier near the detector provides impedance matching so that the pulses can be transmitted along 50 ohm coaxial cable into a 50 ohm terminating impedance. A cable length of up to 10 metres (33 feet) should cause no problems. The cable will of course attenuate the pulse amplitude, with the higher frequency components of the pulse being attenuated more. Assuming that the lowest and dominant frequency component of the pulse is around 20 MHz, corresponding to a typical 25 nsec pulse width, some typical cable attenuations are:

Coax type	Attenuation per 10 m	Distance to reduce amplitude to half
RG-58 (Belden 9201)	0.55 dB	108 m
RG-174 (Belden 9239)	1.18 dB	51 m
RG-316 (M17/172)	1.16 dB	52 m

If you connect the photomultiplier output directly to the C400, then the cable length must be much less. The capacitive load of the cable will tend to integrate the signal, reducing its height and degrading the pulse pair resolution. We recommend maximum capacitive load of 50 pF, at which point the amplitude of the fastest detector pulses is reduced to about 67% when driving a 50 ohm load. This corresponds to a length of about 50 cm of miniature coaxial cable such as RG-174 or RG-316.

The high voltage cable should be coaxial with a DC rating sufficient for the maximum voltage you expect to use. RG-59 is a good choice, with a DC rating of 2.1 kV. However many photomultipliers do not require a voltage this high, and use a smaller coaxial cable such as RG-316 for convenience. SHV connectors that accept this smaller cable size are available from Pasternack (PE4498) and Tyco (5051426-8 and 51426-8).

The pre-amplifier power cables must be well-screened multi-way cable. It is important to avoid injecting noise into wideband pre-amplifiers. We recommend that the screen is connected at the C400 end and the pre-amplifier end, although you may find that it is beneficial to connect at one end only in some noise environments.

Four-channel Pulse Counting Detector Controller

HOW THE C400 WORKS - AN OVERVIEW

The C400 is a very powerful instrument which provides many facilities for measuring fast pulses at high rates on four identical channels. In this section we shall take a general overview of these capabilities.

PULSE HEIGHT DISCRIMINATION AND COUNTING

Incoming fast pulses are terminated in 50 ohms and the resulting voltage is compared to user defined discriminator levels. If the pulse height lies between the low and high discriminator thresholds, then a counter is incremented, and a standard height fast pulse is transmitted from the monitor output. Note that pulses from photomultipliers and other electron multiplier types are negative going. The C400 can also handle positive going pulses if required, to suit inverting pre-amplifiers, or detectors that produce positive pulses.

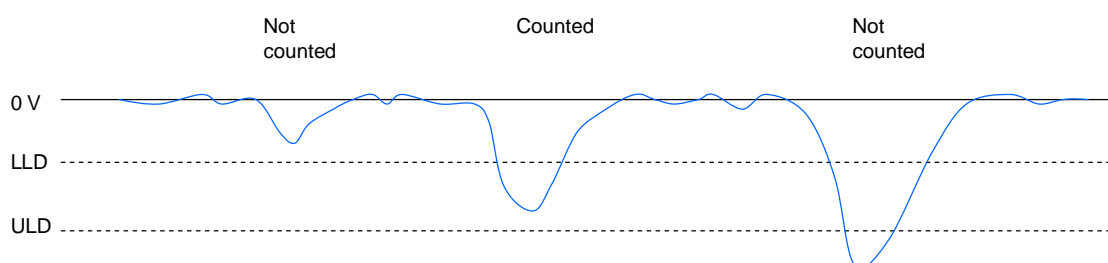


Figure 6. Comparator operation

The comparator circuit includes some hysteresis in the thresholds, to avoid spurious counts due to noise on the trailing side of a pulse.

The pulses are counted for a user-defined integration time, then the counter is reset and the process repeats. There is no deadtime between the integrations. You get better counting statistics for any given integration if the integration time is long, and a reduced data rate, but of course you trade off time resolution of the count rate.

TRIGGERING

In many cases you will need to coordinate the C400 measurements with external events. A TTL signal or on the gate input can trigger a pre-defined acquisition sequence. The gate in signal is echoed quickly on the gate out connector, allowing you to connect multiple units in daisy chain. There are various external triggering modes, from a simple start trigger that starts a predefined sequence of integrations, to modes that allow you to use the gate line to start, stop or pause the counting.

If you don't need to synchronize to external events, you can use internal triggering, which is similar to the auto run facility on an oscilloscope. The C400 generates its own triggers, so counting starts as soon as you initiate an acquisition.

Four-channel Pulse Counting Detector Controller

ACQUISITION MODES

The C400 can generate data very quickly, at a rate higher than can generally be accepted in real time by either the communications link or the host computer. The C400 therefore provides buffering schemes that allow you to capture blocks of data at any rate.

Unbuffered mode does not store any results in the C400. The values are simply transmitted to the host computer as fast as the C400 can send them, and as fast as the host it can accept them. At the shortest integration times, some values may be lost because of limitations in the communications rate, but this is always obvious in any data you take because the readings are individually time-stamped

Buffered mode uses the internal memory of the C400 to store up to 65535 sets of readings from the four channels. They are transmitted to the host in groups of not more than 100, so that the data you see on the host computer is always relatively recent. There is no break in the readings, even at the shortest integration times.

Burst mode subdivides the buffer memory you allocate into smaller bursts of readings. Each burst is accumulated in response to an external trigger. This mode is useful if you have bursts of high speed data, interspersed with periods when there is no useful data.

DETECTOR SERVICES

Each of the four channels provides the services typically needed to operate a detector system. The C400 can be fitted with four individual programmable high-voltage modules to provide bias voltage to the detectors. Most applications will require all four to be specified with the same polarity and voltage range, although the C400 can handle them all being different if necessary. If you are instead using separate high voltage supplies, then you can use the analog output voltages from the C400 to program the voltage from those supplies. ± 12 volt power rails are provided for pre-amplifiers.

SELF-TESTING AND CALIBRATION

Some scintillation detectors include test pulser LEDs. A short pulse of light from the LED causes the photocathode to emit, and you get a signal pulse from the anode of the photomultiplier, which can be discriminated and counted in the normal way. The C400 includes TTL test output pulses for each channel, which you can program with period and pulse width. With a suitable series resistor for the LED, you should count the same pulses coming in as you are putting out, and therefore validate the complete detector and counting system.

The C400 uses precision digital to analog converters (DACs) to control the high voltage supplies and to set the discriminators, and analog to digital converters (ADCs) to readback the high voltages. Calibrations are provided for all these converters so that the high voltage and discriminators are as accurate as possible.

Four-channel Pulse Counting Detector Controller

MONITOR OUTPUTS

Each C400 channel has an associated monitor pulse output. Each valid input pulse that lies between the lower and upper discriminator thresholds is counted and also generates a fixed amplitude output pulse on the corresponding monitor. These pulses may be fed to independent counters, or be used to drive ratemeters or oscilloscope external triggers.

COMMUNICATION TO THE HOST COMPUTER

The C400 provides three alternative communications interfaces to the host computer. You must work with a suitable host computer and software to set up acquisitions on the C400, and read back and display the results.

The most commonly-used is the standard 100/10BaseT Ethernet interface which supports TCP/IP and UDP protocols. The C400 can serve multiple clients on the network.

There is also a 10 Mbps serial fiber optic channel which provides compatibility with the Pyramid Technical Consultants, Inc. fiber optic real-time loop controllers, and thus the complete range of Pyramid products. The fiber optic interface can alternatively be used to connect slave devices to the C400, for example the Pyramid M10 device, to provide extended I/O capability,

Finally, a serial interface which can be operated with RS-232 or RS-485 levels provides a simple interface for applications where high data rates are not needed. An ASCII protocol is supported, based on the SCPI standard. The C400 will support parallel communication on Ethernet and serial interfaces.

The initial C400 software release supports the Ethernet and ASCII serial interfaces. The other interfaces will be added in a future release.

EMBEDDED SOFTWARE

The C400 runs an embedded version of the Linux operating system on NIOS processors implemented in the FPGAs. Four firmware files make up a full release, and you can update with a single file which integrates all the releases, and thus ensures that you have compatible versions.

Four-channel Pulse Counting Detector Controller

GETTING STARTED USING THE PTC DIAGNOSTIC G2 HOST PROGRAM

The PTC DiagnosticG2 is a stand-alone program which allows you to read, graph and log data from the C400, and set all the important acquisition control parameters. It supports communication via any of the interfaces. For some applications it may be adequate for all of your data acquisition needs. The Diagnostic uses the same function library that is exposed for users who develop their own host applications, and therefore also serves as a debugging aid.

PTC DiagnosticG2 was introduced to support the G2 range of Pyramid Technical Consultants, Inc. products, which feature embedded Linux processors and built-in Ethernet interfaces. It is not compatible with the PSI Diagnostic program which supports previous Pyramid products. However the PTC Diagnostic G2 program will be extended in the future to add support for all the previous products.

Your C400 was shipped with a USB memory stick with the installation files you need. We recommend that you copy the files into a directory on your host PC. Check the Pyramid Technical Consultants, Inc. web site at www.ptcusa.com for the latest versions.

INSTALLATION

The program runs under the Microsoft Windows operating system with the 3.5 .NET framework, or later. This has to be installed before the PSI DiagnosticG2. New Windows PCs have .NET already installed. It can be downloaded from the Microsoft web site at no charge.

Copy the installer file PTCDiagnosticSetup-Vx_xx.msi to the hard drive of the host computer, where x_xx is the version of the Diagnostic. The program will run on Windows XP, Vista and Windows 7. The PC must have a standard Ethernet port.

Run the installer and follow the prompts. On Windows 7 systems you will be asked to allow the installer to make changes to the computer after the confirmation stage. You will require administrator access to permit this.

Four-channel Pulse Counting Detector Controller

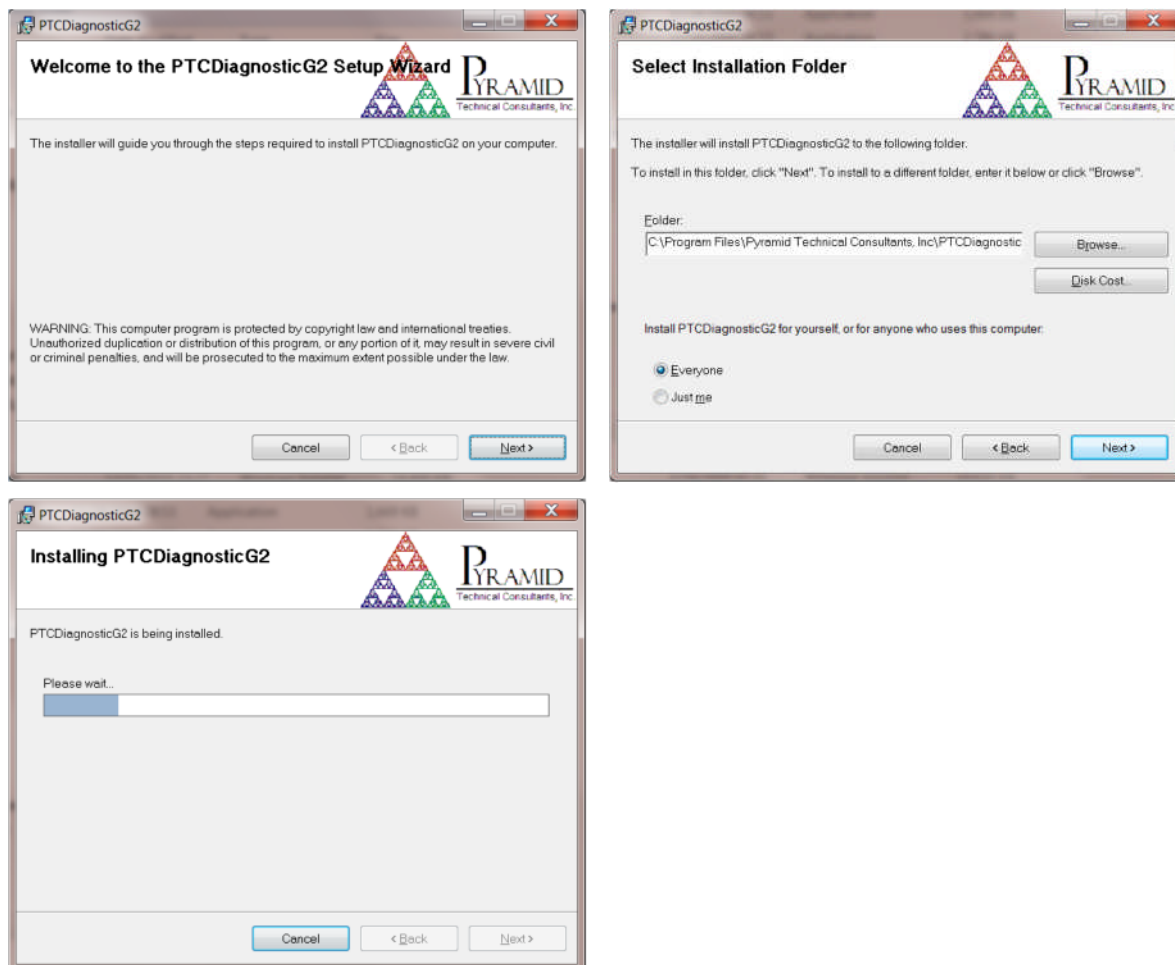


Figure 7. PTC DiagnosticG2 installation

The installer will create a subdirectory in the Program Files directory containing all the executables and configuration files, and create shortcuts on your desktop and in the Start menu.

CONNECTING TO THE C400

The following steps take you through the process of connecting to the device.

- 1) It is simplest to start with a direct connection from your host computer to the C400 using a CAT5 or CAT6 Ethernet cable as shown below. The network cable can be a patch or a crossover type – the C400 automatically adjusts itself to suit.

Four-channel Pulse Counting Detector Controller

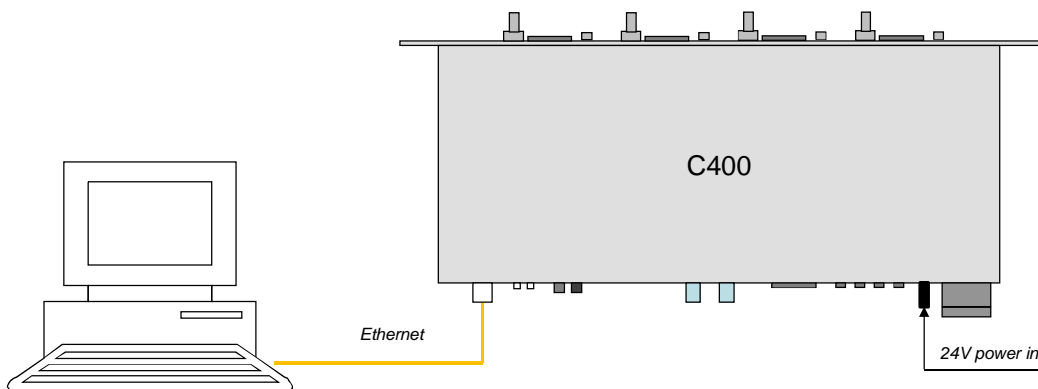


Figure 8. Direct Ethernet connection.

Once you have established reliable communication, and set a suitable unique IP address, then you can move the C400 onto a general local area network if required.

2) The device is set with IP address 192.168.100.20 at shipment. Once you have a connection you can change this setting as required. Set up your host PC Ethernet port with a fixed, non-conflicting IP address in the same subnet range, for example 192.168.100.11.

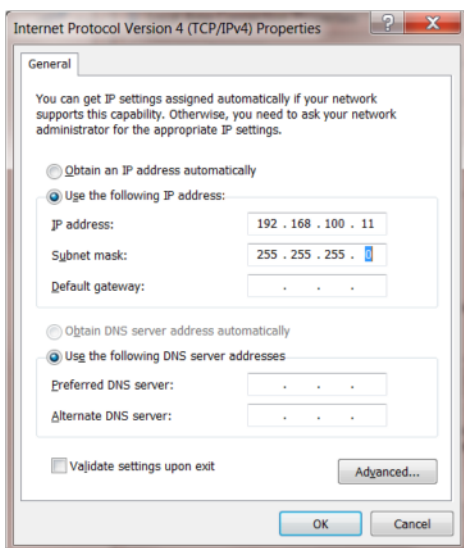


Figure 8a. Configuring the IP address on the host computer.

3) Turn on 24 V DC power to the C400, but make no other connections. The two power LEDs on the rear panel should illuminate when the power is applied, and the cooling fan will start up. While the device is booting, the left hand bank of four LEDs on the rear panel cycle. When the device is ready, the "Active" and power LEDs only should be illuminated.

Four-channel Pulse Counting Detector Controller

4) Make the Ethernet connection from the host PC to the C400. You should see activity on the LEDs that are mounted in the C400 RJ-45 connector. Check that you can ping the device from a command window prompt.

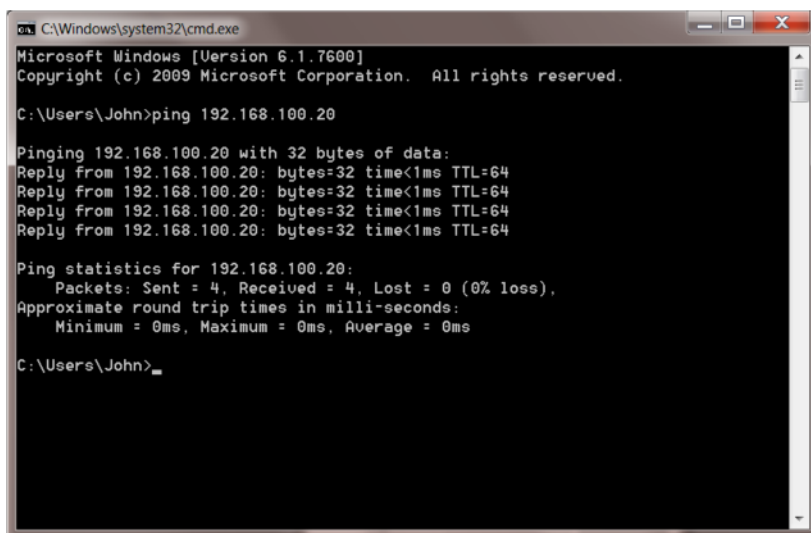


Figure 9. Ping test of the Ethernet connection.

The Windows firewall may block communication with the C400. It is simplest to start by turning off the firewall while you are testing. To maintain security, you can disable any wireless adaptor on the PC while you are doing this. Once you have established communication, then you can try restoring the firewall, and set up permissions for the PTC Diagnostic to communicate through the firewall.

6) Start the PTCDiagnosticG2 software. It will start with the Discover Devices dialog open.

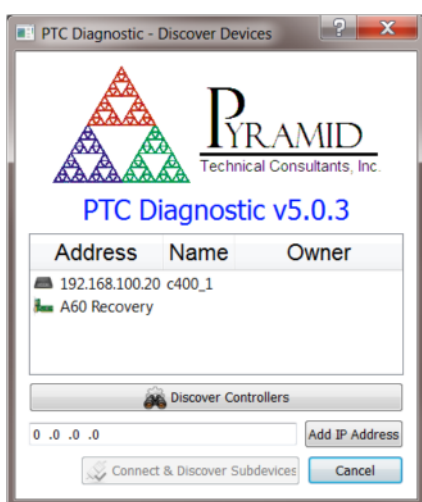


Figure 10. Discover devices

Four-channel Pulse Counting Detector Controller

Click on Discover Controllers and the software will search all IP addresses within its available network looking for compatible Pyramid devices. Since you have the simplest possible network, it should find only the C400 you are working with. In some circumstances the discovered devices may include Pyramid tools such as the A60 recovery utility. You can ignore these.

When you click on the discovered C400, the Connect and Discover Subdevices button is enabled. Click on this to establish the connection to the C400.

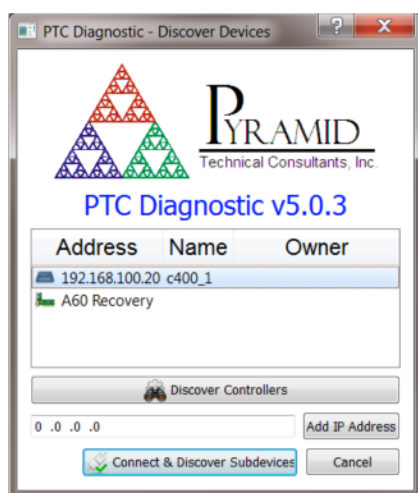


Figure 11. Ready to connect to the C400

A window for the C400 will open, and you will see messages in the message area as the C400 adds your PC as a host.

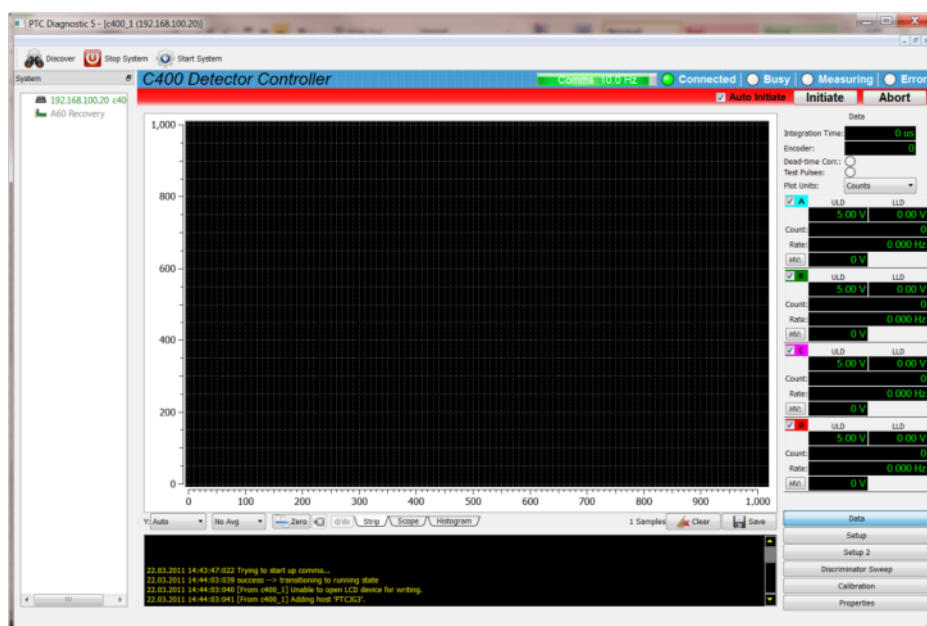


Figure 12. Opening the C400 window

Four-channel Pulse Counting Detector Controller

SCREEN LAYOUT

The C400 user interface screen is divided into two halves plus a top banner area.

The top banner contains the following indicators:

Comms bar	When moving, this indicates that messages from the C400 are being received by the PTCDiagnosticG2. The message frequency is displayed.
Connected LED	When lit, this indicates that communications are valid and the system is not in error.
Busy LED	When lit, this indicates the C400 is busy and cannot respond to inputs, for example while performing a calibration.
Measuring LED	When lit, this indicates that data acquisition has been initiated.
Error LED	When lit, a communication error between the C400 and the host has been detected.

On the left below the banner there is a graphic display of the data with controls for how the data is plotted below it.

Below the graphic is a message window which reports all the commands issued to the C400 by the PTC Diagnostic program, and the corresponding acknowledgements. Generally you can ignore this display, but it will be valuable for diagnosis if you have any operating problems.

On the right there is a screen area which changes according to which display option you select.

Data	A numeric display of the data on each channel, and check boxes to allow you to suppress the graphic display of any subset of the channels. The color codes correspond to the colors of the traces of the strip chart and scope mode graphics. Any channel which is highlighted in red is overrange.
Setup	Controls to set the acquisition mode, deadtime correction, signal polarity, discriminator levels, high voltage settings and data buffering.
Setup 2	Controls to set up test pulse outputs and the position encoder input.
Discriminator Sweep	Controls for the automatic discriminator sweep facility which will create a pulse height spectrum for you, and thus allow you to adjust discriminator levels and high voltage settings to optimum values.
Calibration	Controls for the calibration of all built-in digital to analog and analog to digital converters.
Properties	Firmware version display and update controls, and IP address setting.

Four-channel Pulse Counting Detector Controller

GRAPHIC DISPLAY

There are three ways of displaying incoming data in a graphical way, as a rolling strip chart, as a scope display and as a histogram.

Strip	Data from the four channels, with the selected averaging, is plotted onto a rolling strip chart as it is acquired. The horizontal axis is the time of acquisition and the vertical (y) axis is the current. You can select automatic or fixed vertical scaling.
Scope	A 256 sample data block that has been buffered in the C400 is displayed as a whole. The next available data block replaces the previous one as a whole in the display. This mode is similar to the operation of a digital storage oscilloscope, and is most appropriate for displaying buffered signals acquired in response to an external trigger. The screen looks the same as strip mode.
Histogram	The count or rate in each channel is displayed as a vertical bar. This mode emulates a ratemeter, and can be useful for instrument tuning.

Clicking in the graphic area adds a color-coded cursor which gives the count at that time (strip display) or continuously in that channel (histogram).

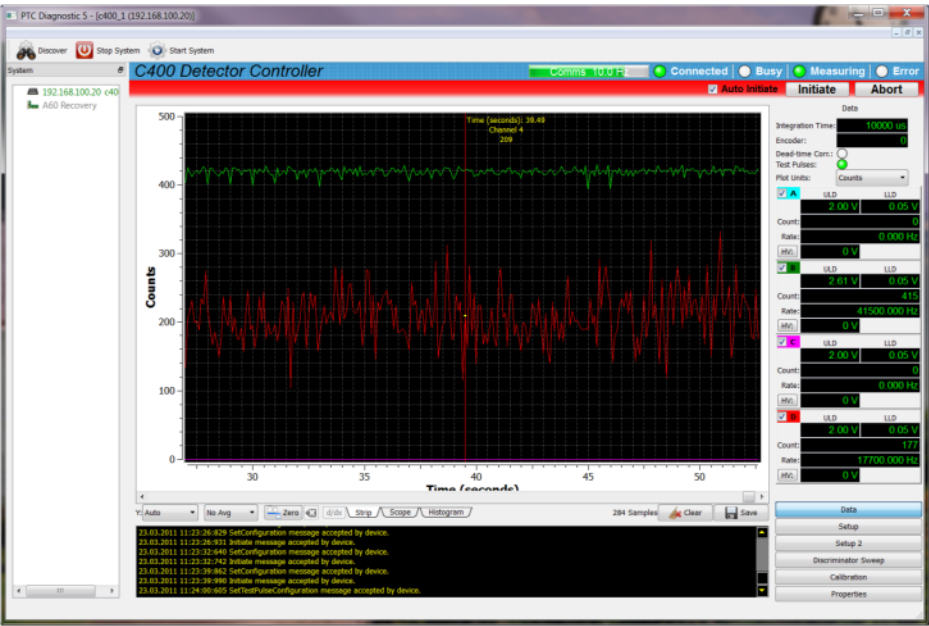


Figure 13. Strip chart data display showing cursor. The scope display looks similar.

Four-channel Pulse Counting Detector Controller

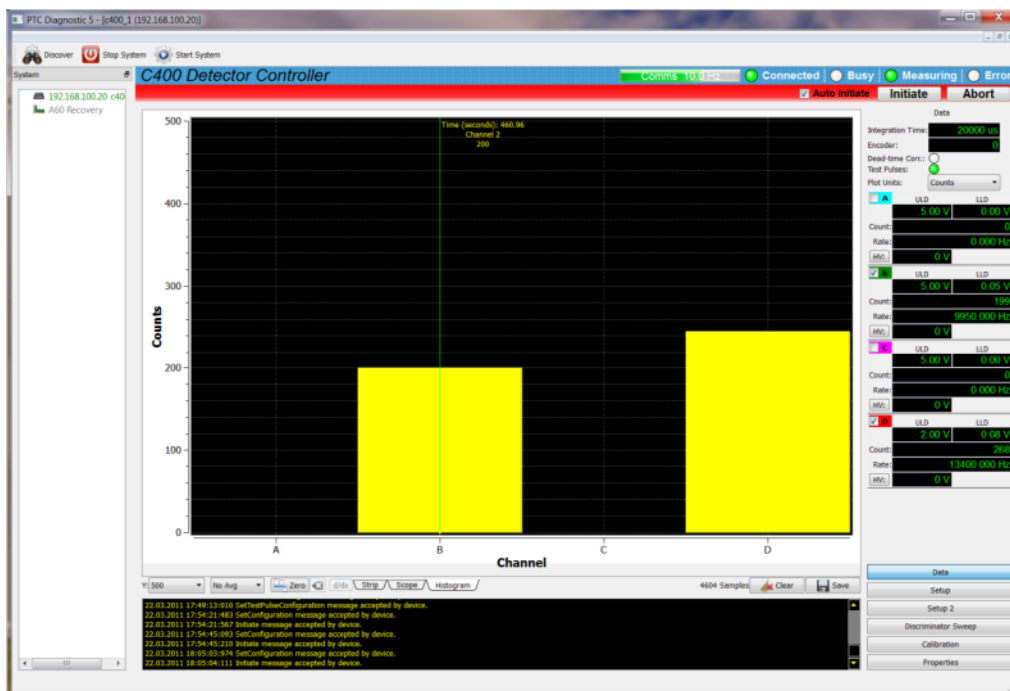


Figure 14. Histogram data display showing cursor.

You can choose whether to display the data as raw counts per integration period, or as the equivalent count rate in Hz or kHz. The displayed data can be filtered by a low pass filter to allow you to pick trends from noisy data.

Four-channel Pulse Counting Detector Controller

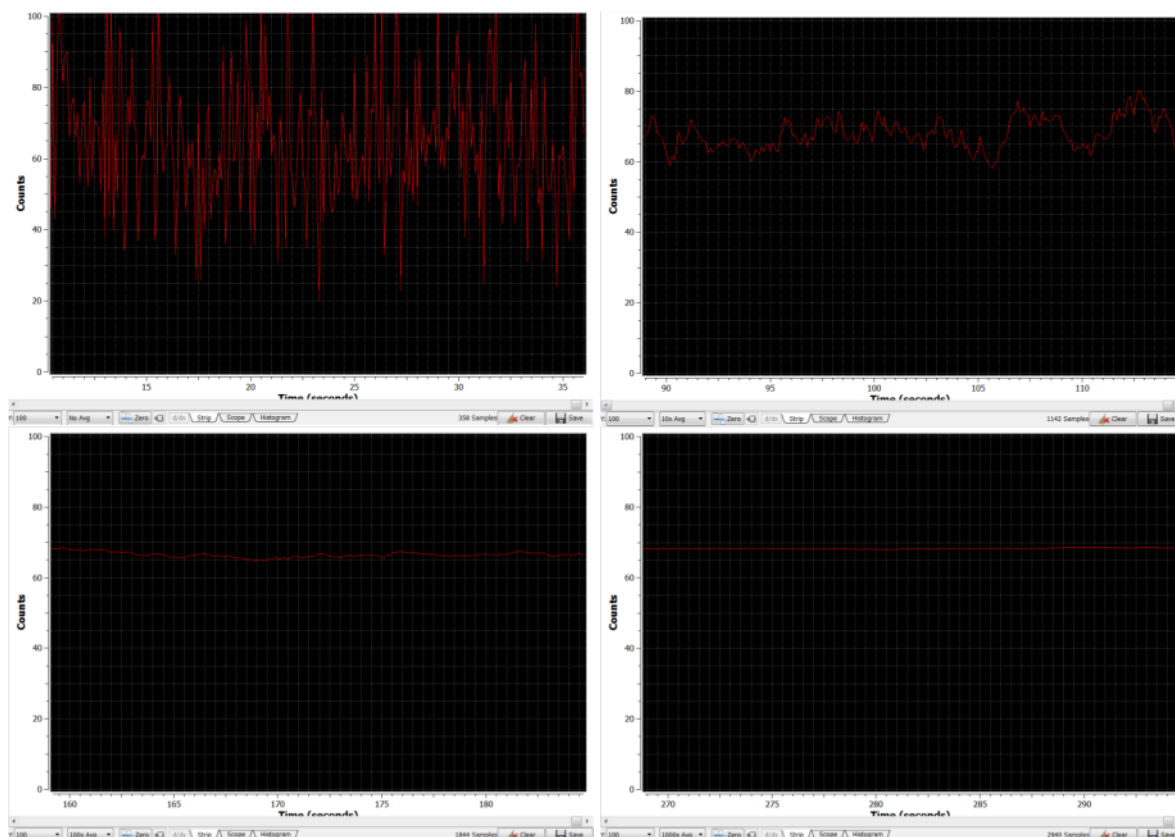


Figure 15. Increasing the low pass filtering of the displayed data.

You can remove any zero offset count rate and thus display changes relative to that offset. This means you may see negative frequencies.



Figure 16. Zero offset subtraction inactive and active.

All of these display options can be selected independently of the acquisition mode, and do not affect the raw data or the logged data, only how it is displayed on the graphic.

Any subset of the four channels can be selected for display on the graphic by checking and unchecking the channels in the Data area. The vertical scale can be set to autoscale, or to be set to various fixed values using the drop down control.

Four-channel Pulse Counting Detector Controller

DATA LOGGING

The PTCDiagnosticG2 software has a data buffer which can accumulate up to 64535 samples, at which point it stops accumulating. Accumulation starts automatically when you click Initiate. You can capture the contents to a .csv format file at any time using the Save button. Pressing the Clear button clears the buffer and restarts the logging.

Note that the logged data is the raw data: any low-pass filtering or zero offset removal does not affect the values you save. However if you have a deadtime correction active, this will be reflected in the data because it is applied in real time to the incoming data. If you need to record data with and without the deadtime correction, and the data can be reproduced, simply toggle the correction off and take a second data set.

DATA DISPLAY

The data display area shows you the numeric values corresponding to the graphic, plus information such as the integration period, the high voltage readings, the discriminator settings in use and whether deadtime correction or the test pulser is active. The high voltage outputs can be enabled and disabled from this display, but if you want to change the values, you must switch to the Setup display.




Figure 17. Data display.

Counts	You can set whether the display of the data is the counts per integration period, or the corresponding count rate in Hz or kHz.
--------	---

Four-channel Pulse Counting Detector Controller

SETUP DISPLAY: DETECTOR, DISCRIMINATOR AND SCALER CONTROL

The primary Setup screen area has controls and parameters for the detectors, the discriminators and the scalers, and thus controls the basic pulse counting function.

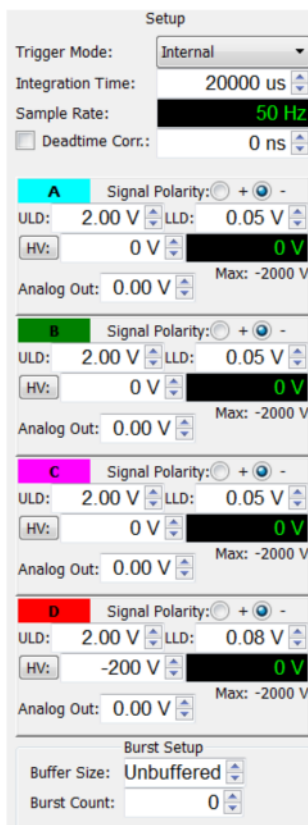

The screenshot displays the 'Setup' screen of the C400. At the top, 'Trigger Mode' is set to 'Internal'. Below it, 'Integration Time' is 20000 us and 'Sample Rate' is 50 Hz. A checkbox for 'Deadtime Corr.' is unchecked, with a value of 0 ns. The screen is divided into four sections for channels A, B, C, and D. Each section includes 'Signal Polarity' (radio buttons for + and -), 'ULD' (Upper Level Discriminator) and 'LLD' (Lower Level Discriminator) thresholds, 'HV' (High Voltage) settings, and 'Analog Out' voltage. Channel A has ULD at 2.00 V, LLD at 0.05 V, HV at 0 V, and Analog Out at 0.00 V. Channel B has ULD at 2.00 V, LLD at 0.05 V, HV at 0 V, and Analog Out at 0.00 V. Channel C has ULD at 2.00 V, LLD at 0.05 V, HV at 0 V, and Analog Out at 0.00 V. Channel D has ULD at 2.00 V, LLD at 0.08 V, HV at -200 V, and Analog Out at 0.00 V. All channels have a 'Max: -2000 V' limit. At the bottom, the 'Burst Setup' section shows 'Buffer Size' as 'Unbuffered' and 'Burst Count' as 0.

Figure 18. Primary setup parameters.

In order for any change to be sent to the C400, you must either press the Enter key on your computer, or click into another parameter field. Any change to the acquisition mode, or to any acquisition parameters, will abort the acquisition in progress. If you have the Auto Initiate box checked, a new acquisition with the new parameters will start immediately.

The setup parameters you enter can all be saved to non-volatile memory in the C400 and will be reloaded if you reboot the device.

Four-channel Pulse Counting Detector Controller


Trigger mode	<p>Several trigger modes are available with this drop-down control.</p> <p><i>Internal:</i> The C400 generates its own triggers, so acquisition starts as soon as the unit is initiated</p> <p><i>External Start:</i> After initiation, the C400 waits for a rising (falling) edge on the gate input to start integrating counts</p> <p><i>External Start/Stop:</i> After initiation, the C400 waits for a rising (falling) edge on the gate input to start integrating counts. A falling (rising) edge stops the acquisition.</p> <p><i>External Start/Hold:</i> After initiation, the C400 waits for a rising (falling) edge on the gate input to start integrating counts. A falling (rising) edge suspends the acquisition until the next rising edge.</p> <p><i>External Windowed:</i></p> <p>The way that acquisitions are triggered, and the numbers of readings that are taken in response to triggers is controlled by the trigger mode, the buffer size and the burst count. This is explained in more detail in section 12. If you simply want to acquire counts right away without worrying about the complexities of triggering, select Internal trigger mode and 0 buffer size (unbuffered).</p>
Integration time	<p>This sets the time intervals over which the scalers accumulate counts, and thus the rate that new readings are taken. You can set from 10 μsec to 999999 μsec, You would generally use longer times for low count rates, but note that you can always sum any number of consecutive readings in host-computer data processing to achieve the effect of very long integrations.</p>
Deadtime correction	<p>The C400 can apply a non-paralyzable deadtime correction in real time to the incoming counts in all channels. You can set the deadtime duration in nsec and independently toggle the correction on or off. See section 15 for further details.</p>
Signal polarity	<p>The C400 can handle negative-going or positive-going pulses, and you can set this polarity independently for each channel.</p>
ULD and LLD	<p>The upper and lower levels for the window discriminator can be set independently for each channel. The polarity of the levels is set automatically according to the signal polarity selection.</p>
HV	<p>If you have the high voltage options installed, you can set the output voltages individually for each channel. The "Max: " value in the display shows the highest voltage and polarity available with the installed option, which could be different for each channel. You must enter a value in range, and with the</p>

Four-channel Pulse Counting Detector Controller

	explicit polarity (that is, if you have negative modules installed, you must enter "-525" to get a negative 525 volt output. The high voltage is only enabled when you press the HV: button; you should then see the voltage ramp to the set point.
Analog out	The C400 has a precision +/- 10 V analog output for each channel. You can set the value directly in this edit box. This may be used for example to program a remote high voltage supply.
Buffer size	If you require contiguous readings at high rates, then you must buffer the data in the C400 because the communication link and the host computer cannot be guaranteed to keep up with the data rate. You can define a buffer of up to 65536 readings. Data is uploaded to the host in blocks of 100 while the acquisition is in progress. If you set a buffer of 0, then buffering is disabled and the data is sent to the host computer as fast as it can be accepted.
Burst count	The burst count is the number of readings that are taken after each trigger. The buffer may be considered as a reserve from which fixed portion is drawn whenever there is a trigger.

SETUP 2 DISPLAY: TEST PULSER AND ENCODER INPUT CONTROLS

The second setup page lets you configure the test pulser outputs and the encoder input.

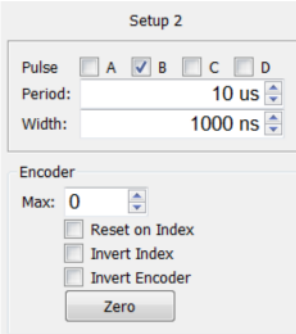

The screenshot shows a software window titled "Setup 2". It contains two main sections: "Pulse" and "Encoder".

Pulse section:

- Four checkboxes labeled A, B, C, and D. Checkboxes A, B, and C are unchecked, while checkbox D is checked.
- A "Period:" label followed by a text box containing "10 us" and a small up/down arrow icon.
- A "Width:" label followed by a text box containing "1000 ns" and a small up/down arrow icon.

Encoder section:

- A "Max:" label followed by a text box containing "0" and a small up/down arrow icon.
- Three checkboxes: "Reset on Index", "Invert Index", and "Invert Encoder". All three are unchecked.
- A "Zero" button located below the checkboxes.


Figure 19. Secondary setup parameters.

Four-channel Pulse Counting Detector Controller

Pulse setup	You can enable the diagnostic pulse output individually for the four channels by checking any combination of the boxes A to D. This facility is typically used in scintillation detectors that include a test LED in the housing, to confirm the complete detector and counting system is working. The pulse width is adjustable down to 1 nsec and the period down to 1 μ sec (corresponding to 1 MHz).
Encoder	You can limit the encoder count to a specific maximum, beyond which it starts again from zero. This is typical for shaft encoders. Setting zero means there is no limit to the count. If there is an index signal available, you can cause this to zero the counter, and you can choose the polarity of the index. The counting sense relative to the physical motion sense can be swapped by checking the Invert Encoder box.

DISCRIMINATOR SWEEP

In order to set the detector high voltage and the discriminator levels, you need information about the distribution of pulse heights. The C400 provides an automatic discriminator sweep facility that displays a spectrum of pulse heights, allowing you to make informed decisions about the settings. The discriminator and high voltage settings are reproduced in this display area so that you can enter the values you choose conveniently. See section 14 for further details.

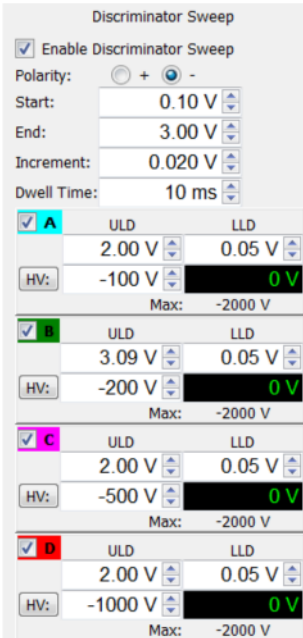

Discriminator Sweep

☒ Enable Discriminator Sweep

Polarity: ☐ + ☒ -

Start: 0.10 V

End: 3.00 V

Increment: 0.020 V

Dwell Time: 10 ms

	ULD	LLD
<input checked="" type="checkbox"/> A	2.00 V	0.05 V
HV:	-100 V	0 V
Max:	-2000 V	
<input checked="" type="checkbox"/> B	3.09 V	0.05 V
HV:	-200 V	0 V
Max:	-2000 V	
<input checked="" type="checkbox"/> C	2.00 V	0.05 V
HV:	-500 V	0 V
Max:	-2000 V	
<input checked="" type="checkbox"/> D	2.00 V	0.05 V
HV:	-1000 V	0 V
Max:	-2000 V	

Figure 20. Discriminator sweep controls.

Four-channel Pulse Counting Detector Controller

Enable discriminator sweep	Discriminator sweep is a special acquisition mode that you must enable by checking the box. The usual upper and lower levels are overridden by the dynamic values used in the sweep.
Start, End	These controls set the starting and ending lower discriminator levels for the sweep.
Increment	The window that is swept through the pulse height range, and thus the resolution of the measurement, is set by the increment control. The steps are the same size, so that the resulting spectrum is contiguous.
Dwell time	This sets how long data is accumulated at each step.

CALIBRATION

The C400 stores calibrations for all its DACs and ADCs to improve their absolute accuracy. Any of the values can be edited on this display, and are stored in non-volatile memory when the values are entered. You should clearly be careful not to enter incorrect values as you will get confusing results. The discriminator upper and lower level voltages and the high voltage readback can be calibrated automatically, either selectively by channel, or altogether.

automatically, either selectively by channel, or altogether.

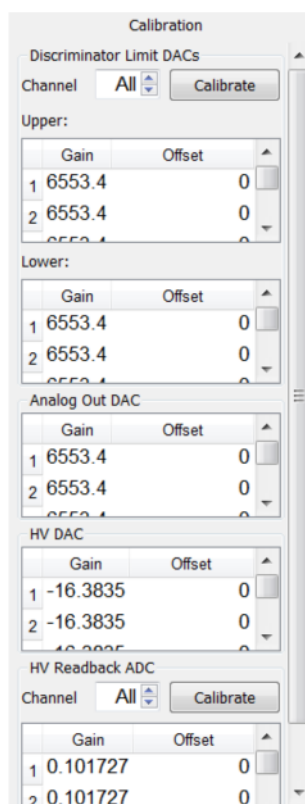


Figure 21. Calibration control and display.

Discriminator	The gain values are in DAC bits per volt, and the offsets are in volts. The
----------------------	---

Four-channel Pulse Counting Detector Controller

limit DACs	Calibrate button at the top of the display adjusts the offsets automatically. The offset can be adjusted to include the offset voltage of any pre-amplifier that is attached to the channel.
Analog out DAC	The gain values are in DAC bits per volt, and the offsets are in volts. These values can be altered manually based upon measurement of the analog output voltage with an independent meter.
HV DAC	The gain values are in DAC bits per volt, and the offsets are in volts. These values can be altered manually based upon measurement of the high voltage output with an independent meter.
HV readback ADC	The gain values are in volts per ADC bit, and the offsets are in volts. These values are adjusted automatically by the C400 when you press the Calibrate button at the bottom of the display area. You should ensure that there is nothing connected to the high voltage connectors before pressing Calibrate. The process takes a few minute to complete; you will see the Busy indicator lit. After calibration you should see good agreement between the commanded HV values and the readback.
Clear All Calibrations	Clicking the Clear button deletes all the stored calibration values, setting all the offsets to zero and all the gains to default values.

PROPERTIES

Firmware update

The C400 has four embedded firmware installations, which need to be compatible with each other. The installations are the Linux operating system, the NIOS processor real time application and the FPGA programs for primary and secondary FPGAs. In order to simplify the task of updating the firmware and ensuring the versions are all compatible, Pyramid Technical Consultants, Inc. releases C400 firmware as a single integrated file which is unpacked and uploaded to the correct memory areas by the software.

You should ensure that the C400 and the host computer will not be disturbed for the four minutes approximately that the update requires. To perform the update, click the Update All Firmware button on the Properties screen area. The software will now warn about the update process, and ask if you want to proceed.

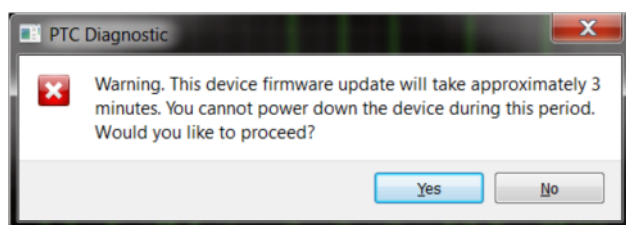


Figure 22. Firmware update warning.

Four-channel Pulse Counting Detector Controller

Navigate to the update file, which is typically about 5 Mb in size and has filename C400_Firmware_x.x.x.x.zip where x.x.x.x is the version number.

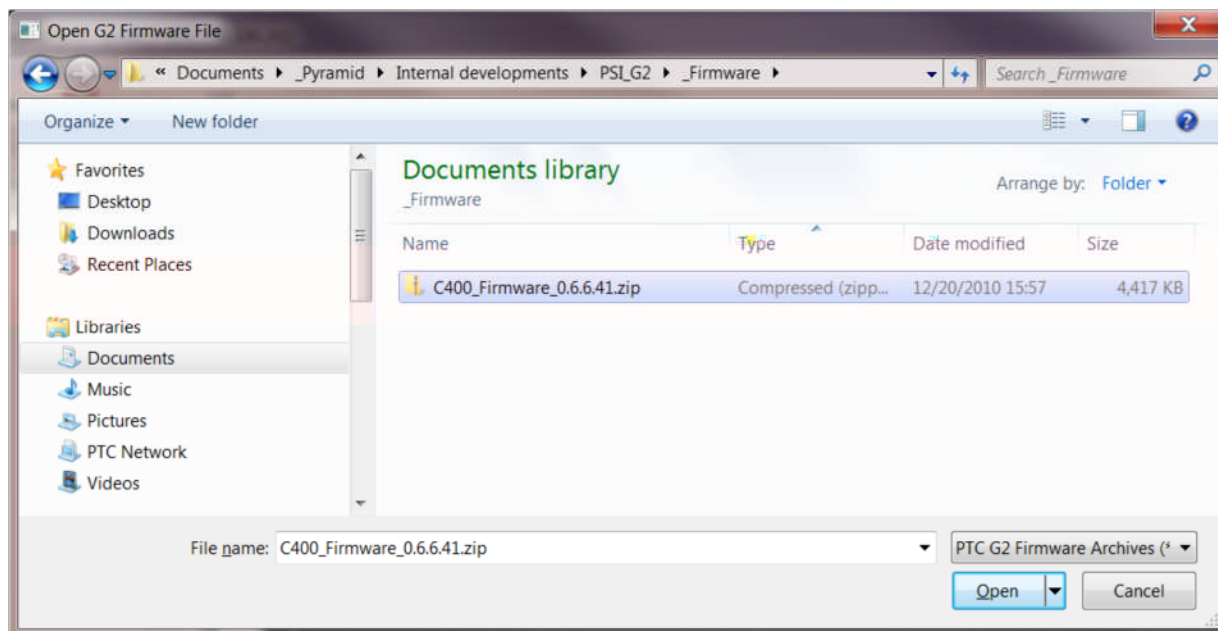


Figure 23. Selecting the firmware file.

While the update is taking place, you will see a countdown of the time remaining, and you can see detailed activities reported in the message window.

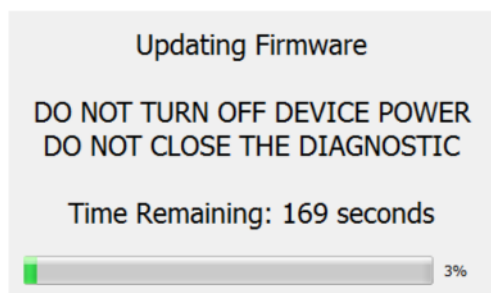


Figure 24. Firmware update status window.

If the process does get interrupted, it is possible that the C400 will not be able to do another update because core programs have become corrupted. In this case you should contact Pyramid or your supplier who can organize a recovery. When the update is complete, you should reboot the device by power cycling it, and check the Properties screen area to ensure that the new versions are correctly reported.

Four-channel Pulse Counting Detector Controller

Changing the IP address

You can alter the C400 fixed IP address to any valid setting.

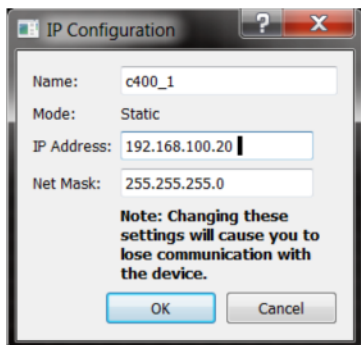
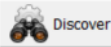


Figure 25. IP address configuration dialog.

Note that the address change takes effect immediately, so you will lose communication in the session you were running. Simply click  on the Diagnostic to rediscover the device. If you have set the IP address outside the subdomain that is accessible to your host computer, you will need to change the IP address of the computer so that it is in the same subdomain as the C400.

Restoring the default IP address

You can force the C400 to restore its default IP address of 192.168.100.20 by power cycling and keeping the rear panel reset button pressed while it boots up.

Four-channel Pulse Counting Detector Controller

C400 CIRCUIT OVERVIEW

DISCRIMINATOR INPUTS

The C400 uses fast level discriminators to compare the size of the signal at its input with thresholds set by digital to analog converters. The signal is terminated in 50 ohms and fed in parallel to the lower and upper level discriminators. The logic outputs are delivered to a high-performance FPGA (Field Programmable Gate Array).

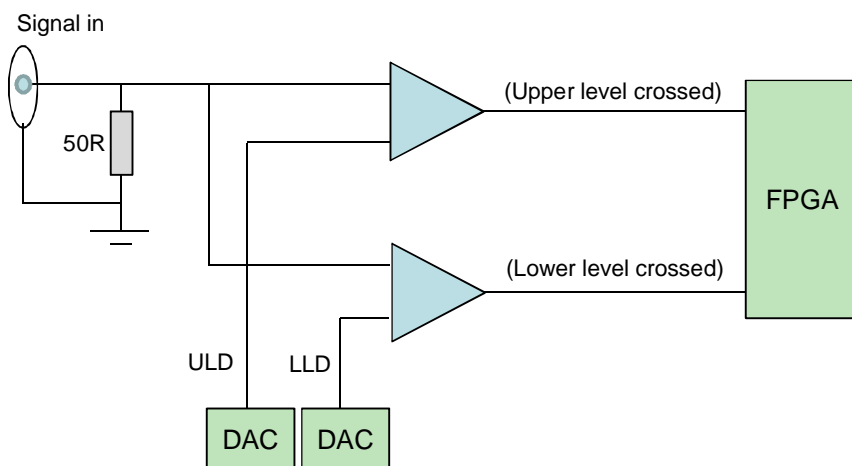


Figure 26. The basic discriminator circuit.

The discriminators have hysteresis of about 30 mV to prevent multiple triggering on noise. The C400 sets the levels as necessary to handle negative or positive going pulses, as selected by the user.

The FPGA sets a bit when the low level discriminator produces an output. If the upper level also produces an output, then the bit is reset, and nothing is counted. If the bit is not reset, it is counted into the scalers. The scalers have a depth of 32 bits, or 4.3 billion counts. They are read out transparently, with no interruption to the counting.

OVERALL CIRCUIT

The C400 is based around high-performance two field programmable gate arrays (FPGAs) which handle all data collection, processing and communications. One is on the standard Pyramid A60 processor module, which includes the NIOS II embedded processor core which runs the C400 real time application. The A60 manages all the communications with the host computer. It reads the communication mode switches, jumpers on the main board that inform it about system configuration, and drives the system status LEDs.

Four-channel Pulse Counting Detector Controller

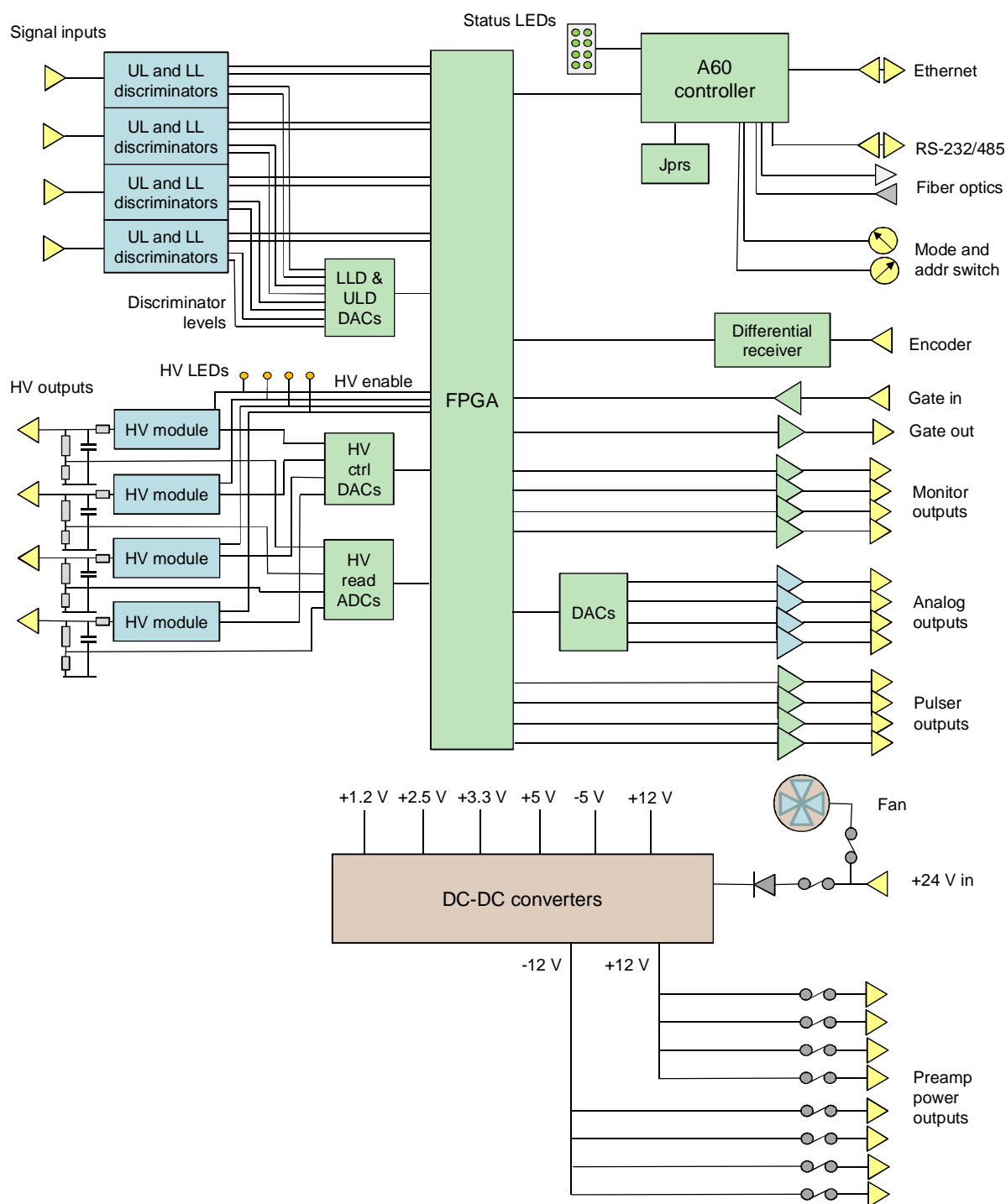


Figure 27. C400 block schematic.

The second FPGA is on the main board, and handles all the specific tasks such as discriminator logic, scalers, encoder pulse counting, and all DAC and ADC management.

Four-channel Pulse Counting Detector Controller

The high voltage outputs include 100 μ sec RC filters to reduce ripple. The high voltage readback is taken from resistive dividers directly on the outputs.

24 VDC input power enters via a 1.1 A resettable fuse, and is used directly to power the cooling fan. The 24 V input is protected against polarity reversal by a series diode, and from transients by series inductors and transorbs to chassis which limit excursions to 6 V (not shown on the block schematic). DC-DC converters generate the ± 12 VDC supplies for the pre-amplifier outputs, +12 VDC for the high voltage modules, ± 5 VDC for the discriminators, and +3.3 VDC, +2.5 VDC and +1.2 VDC for the digital circuits. The pre-amplifier power outputs are individually protected by 200 mA resettable fuses.

Four-channel Pulse Counting Detector Controller

HIGH VOLTAGE SUPPLIES

SETTING THE HIGH VOLTAGE SUPPLIES

The C400 is available with a one-watt high voltage supply suitable for biasing photomultipliers, continuous dynode multiplier and diodes. The voltage range can be specified at time of purchase from 200, 500, 1000 and 2000V with either polarity. This can be done individually by channel, although this is not recommended due to the risk of connecting a detector to the wrong high voltage output. The A60 processor reads internal jumpers to detect the supply rating. Note that the supplies are not intended to deliver output voltage less than about 10% of their maximum rating.

The front panel HV on LED illuminates when the adjacent output is enabled. The set value can be adjusted at any time. If you are using the PTC Diagnostic software and have checked Auto Initiate, any measurement in progress will restart so you can see the effect of the change at once.

The outputs of the high voltage modules are filtered by an RC filter to reduce ripple and noise. For large current drains there will be a voltage drop across the series resistor, but you can see the actual output voltage because readback is taken from a voltage divider after the filter, directly on the HV output. The voltage divider places a fixed 60 Mohm load on the supply.

The following table gives the smallest value of external load resistance that the various supply options can drive at full voltage, the filter time constant and also the volts dropped across the series resistor at maximum current.

<i>HV module</i>	<i>Current rating</i>	<i>Smallest load resistance</i>	<i>Filter time constant</i>	<i>Filter series resistor</i>	<i>Voltage drop across filter no load</i>	<i>Voltage drop across filter full load</i>
+2000, -2000	0.5 mA	4.3 Mohm	110 μ sec	33.2 kohm	1 V	17 V
+1000, -1000	1.0 mA	1.0 Mohm	33 μ sec	10 kohm	0.2 V	10 V
+500, -500	2.0 mA	250 kohm	16 μ sec	4.7 kohm	< 0.1 V	10 V
+200, -200	5.0 mA	40 kohm	0	0	0	0

If the readback value differs from the setpoint by more than the expected drop across the filter resistor, you know that the output is either being overloaded by a low resistance to ground, or that it is being driven by another source of higher compliance. Be aware of how the high voltage readback calibration might affect this. We recommend that the calibration is done with no connection to the HV outputs. Thus only the no-load voltage drop across the filter is hidden by the calibration, and when you connect a load you should expect to see small reductions in the measured voltage, up to the maxima given in the table. The measured voltage is the actual voltage at the output connector.

Four-channel Pulse Counting Detector Controller

The C400 is available with higher power high voltage modules to special order. Please contact Pyramid Technical Consultants or your supplier for more details.

CAUTION



Do not connect external power supplies to the C400 external high voltage outputs that will drive the built-in supplies away from the voltages they are trying to regulate, or you may cause damage to the C400.

CHANGING THE HIGH VOLTAGE SUPPLY RANGE AND POLARITY

The range and polarity of the high voltage supplies is fixed and must be specified at time of purchase. Units may be returned to the factory to change the high voltage modules if necessary. It not recommended that users change the high voltage supply modules. The jumper settings are given here for reference only. There is a jumper bank for each channel.

























 + HV0  HV1  HV2	-2000 V	 + HV0  HV1  HV2	+2000 V
 + HV0  HV1  HV2	-1000 V	 + HV0  HV1  HV2	+1000 V
 + HV0  HV1  HV2	-500 V	 + HV0  HV1  HV2	+500 V
 + HV0  HV1  HV2	-200 V	 + HV0  HV1  HV2	+200 V

Figure 28. High voltage jumper settings

All other settings are interpreted as no HV option fitted.

Four-channel Pulse Counting Detector Controller

TRIGGERS AND DATA BUFFERING

C400 ACQUISITION STATES

The best way to understand the ways that the C400 trigger modes and buffering options is to consider the internal states of the device concerning data acquisition, and the way you can control the transitions between these states by selecting the trigger modes, the buffer size and the burst count.

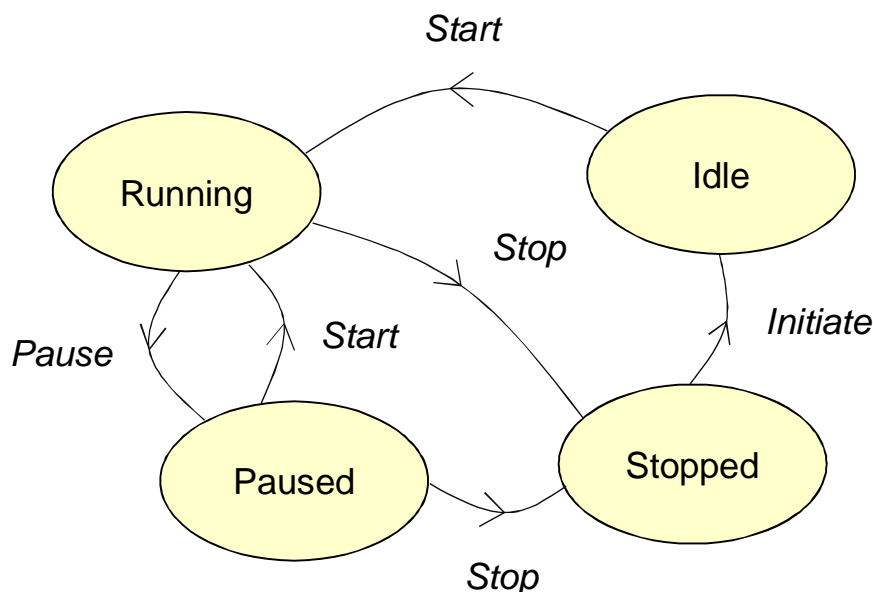


Figure 29. C400 acquisition state diagram.

The C400 is Stopped when there is no data acquisition in progress or pending. When you send the Initiate command it moves to the Idle state. It will then transition to the Running (acquiring data) state when the correct conditions are met. From Running the C400 may then transition to a Paused state while it waits for a condition to be met before resuming, or to the Stopped state to complete the acquisition.

There are three sources regulating the transitions, namely the Start source, the Burst (or Pause) source, and the Stop source. Each can be set to Internal (a C400 internal condition is met) or BNC (an external trigger signal is detected at the gate in connector).

The C400 will transition from Idle to Running if:

- The start source is set to BNC, and the appropriate edge is seen on the BNC gate
- If start source is set to internal, it will automatically go to Running when initiated

Four-channel Pulse Counting Detector Controller

The C400 will transition from Running to Paused if:

- Burst Count is greater than zero and that number of samples is acquired
- The Burst Source is set to BNC, and the appropriate edge is seen on the BNC gate

The C400 will transition from Running or Paused to Stopped if:

- Buffer Size is not Unbuffered (ie, it is greater than 0) and that number of samples is acquired
- An abort command is received
- The stop source is set to BNC, and the appropriate edge is seen on the BNC gate

TRIGGER MODES

In the Custom trigger mode, you have full control over whether state transitions are driven by internal or external sources, for starting, pausing and stopping the acquisition. In the other trigger modes, some of the sources for various transitions are preset, as described in the following table. In the table, NBuf is the buffer size and NBst is the burst size. For every trigger mode you can force the C400 to the stopped state at any time by sending the Abort command.

Mode	Start	Pause	Stop	Notes
Internal	Internal	n/a	n/a	Acquisition will start immediately you send Initiate, and continue indefinitely if unbuffered, or to the lesser of NBuf size and NBst if buffered.
External Start	BNC	n/a	n/a	Acquisition will start when a valid trigger edge is seen after you send Initiate, and continue indefinitely if unbuffered. If buffered, it will do NBst readings on each valid trigger edge until it has acquired NBuf readings.
External StartStop	BNC	n/a	BNC	Acquisition will start when a valid trigger edge is seen after you send Initiate. If unbuffered, readings will continue until the opposite polarity trigger edge is seen. If buffered, it will do , it will stop on the lesser of NBst or Nbuf readings, or on the opposite polarity trigger edge, whichever comes first.
External StartHold	BNC	n/a	n/a	NBst is forced to a value of 1 in this mode. A single reading is taken for each valid trigger edge. This will continue indefinitely if unbuffered. If buffered, it will continue until NBuf is reached.

Four-channel Pulse Counting Detector Controller

External Windowed	BNC	BNC	n/a	Acquisition will start when a valid trigger edge is seen after you send Initiate. They continue until either NBst is reached, or the opposite polarity edge is seen, at which point the acquisitions pause. They resume when the next trigger edge is seen. In unbuffered mode, this continues indefinitely. In buffered mode, it continues until NBuf is reached.
Discriminator Sweep	Internal	n/a	n/a	Special internal mode used for discriminator sweeps only.

BUFFER MEMORY

If you are taking measurements at a relatively low data rate, or if you don't care about having time contiguous data, then you can use the unbuffered mode and use the memory of your PC to accumulate the results. You can tell the maximum data rate before you will start to lose samples if you are using the PTC Diagnostic host software by observing the Comms rate indication. For example, if the rate is shown as 10 Hz, then you can take readings with integration time down to 100,000 μ sec before you will start missing data points in the logged data. The fact that data is missing is evident from the gaps in the trigger counts, as shown in the following examples which used internal trigger mode.

timestamp	triggercount	channel_1	channel_1	channel_2	channel_2	channel_3	channel_3	channel_4	channel_4	encodercount
25.6	0	0	-0.05	0	-0.05	0	-0.05	4357	-0.05	20
25.7	1	0	-0.05	0	-0.05	0	-0.05	4147	-0.05	20
25.8	2	0	-0.05	0	-0.05	0	-0.05	4431	-0.05	20
25.9	3	0	-0.05	0	-0.05	0	-0.05	4124	-0.05	20
26	4	0	-0.05	0	-0.05	0	-0.05	4661	-0.05	20
26.1	5	0	-0.05	0	-0.05	0	-0.05	4559	-0.05	20
26.2	6	0	-0.05	0	-0.05	0	-0.05	5180	-0.05	20
26.3	7	0	-0.05	0	-0.05	0	-0.05	4200	-0.05	20
26.4	8	0	-0.05	0	-0.05	0	-0.05	4405	-0.05	20
26.5	9	0	-0.05	0	-0.05	0	-0.05	5103	-0.05	20
26.6	10	0	-0.05	0	-0.05	0	-0.05	5080	-0.05	20
26.7	11	0	-0.05	0	-0.05	0	-0.05	4433	-0.05	20
26.8	12	0	-0.05	0	-0.05	0	-0.05	4120	-0.05	20

Figure 30. Unbuffered logged data with a 10 Hz comms rate, 100000 μ sec integration time.

timestamp	triggercount	channel_1	channel_1	channel_2	channel_2	channel_3	channel_3	channel_4	channel_4	encodercount
0.01	1	0	-0.05	0	-0.05	0	-0.05	405	-0.05	20
0.11	11	0	-0.05	0	-0.05	0	-0.05	321	-0.05	20
0.21	21	0	-0.05	0	-0.05	0	-0.05	351	-0.05	20
0.31	31	0	-0.05	0	-0.05	0	-0.05	360	-0.05	20
0.41	41	0	-0.05	0	-0.05	0	-0.05	425	-0.05	20
0.51	51	0	-0.05	0	-0.05	0	-0.05	495	-0.05	20
0.61	61	0	-0.05	0	-0.05	0	-0.05	541	-0.05	20
0.71	71	0	-0.05	0	-0.05	0	-0.05	327	-0.05	20
0.81	81	0	-0.05	0	-0.05	0	-0.05	382	-0.05	20
0.91	91	0	-0.05	0	-0.05	0	-0.05	589	-0.05	20
1.01	101	0	-0.05	0	-0.05	0	-0.05	477	-0.05	20

Figure 31. Unbuffered logged data with a 10 Hz comms rate, 10000 μ sec integration time.

Four-channel Pulse Counting Detector Controller

By contrast, the C400 internal data memory can accept data from the counters at the maximum rate, so you should use buffering in cases where it is critical to take time contiguous data.

timestamp	triggercou	channel_1	channel_1	channel_2	channel_2	channel_3	channel_3	channel_4	channel_4	encodercount
0	0	0	-0.05	0	-0.05	0	-0.05	1	-0.05	20
1.00E-05	1	0	-0.05	0	-0.05	0	-0.05	0	-0.05	20
2.00E-05	2	0	-0.05	0	-0.05	0	-0.05	0	-0.05	20
3.00E-05	3	0	-0.05	0	-0.05	0	-0.05	0	-0.05	20
4.00E-05	4	0	-0.05	0	-0.05	0	-0.05	1	-0.05	20
5.00E-05	5	0	-0.05	0	-0.05	0	-0.05	0	-0.05	20
6.00E-05	6	0	-0.05	0	-0.05	0	-0.05	0	-0.05	20
7.00E-05	7	0	-0.05	0	-0.05	0	-0.05	0	-0.05	20
8.00E-05	8	0	-0.05	0	-0.05	0	-0.05	0	-0.05	20
9.00E-05	9	0	-0.05	0	-0.05	0	-0.05	0	-0.05	20
0.0001	10	0	-0.05	0	-0.05	0	-0.05	0	-0.05	20
0.00011	11	0	-0.05	0	-0.05	0	-0.05	0	-0.05	20

Figure 32. Buffered logged data with a 10 Hz comms rate, 10 μ sec integration time.

Four-channel Pulse Counting Detector Controller

MONITOR OUTPUTS

The monitor outputs provide a standard 10 nsec TTL logic pulse into a 50 ohm load whenever an incoming pulse into the corresponding channel falls into the discriminator window. The lag between the arrival of the pulse and the output of the monitor is around 80 nsec after the discriminator conditions are met. In the following figure, the input pulse and the discriminator levels are shown in yellow, and the resulting monitor output pulse in blue.

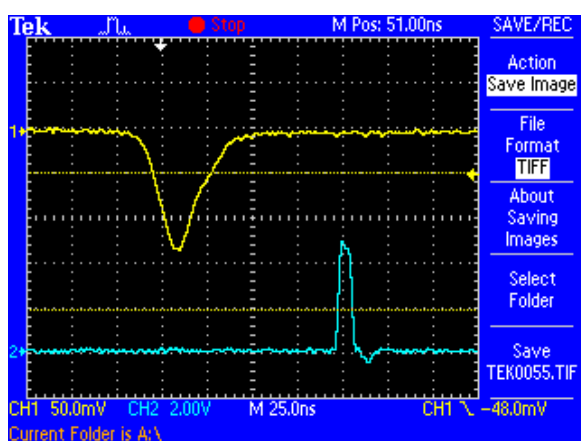


Figure 33. Example of monitor output pulse.

You can use the monitor outputs to drive external scalars, or diagnostic devices such as ratemeters.

If you are viewing the input pulses on a digital oscilloscope, and also feeding them to the C400, then the monitor can make a very useful trigger that selects only the pulses in your range of interest for display. This is not possible on an analog oscilloscope because the monitor occurs after the incoming pulse.

Four-channel Pulse Counting Detector Controller

DISCRIMINATOR SWEEP

In order to make informed decisions about setting detector high voltage and the discriminator levels, you need to know about the range of pulse heights coming from the detector. The C400 provides an automated discriminator sweep function. A small discriminator window is stepped from a starting voltage to an ending voltage that you nominate, and the counts measured at each step are plotted as a function of the lower limit of the window. The sweep repeats so that the counting statistics improve. Making the increment small increases the resolution of the pulse height analysis, but slows down the acquisition. If your count rate is relatively low, it is beneficial to use buffered mode to maximize the number of pulses that you count in a given time— simply set the buffer to the maximum size of 65536.

Although the C400 is not primarily designed to be a pulse height analyzer, you can obtain recognizable energy spectra if the incoming pulses are sufficiently well-defined. This means that you should use a pre-amplifier with sufficient integration to give smooth pulses. See section 16 for more details. The following example shows a spectrum from a Na-22 source, as measured with a LaBr3 scintillation detector, a CP10-A preamplifier and the C400. The 511 keV positron annihilation peak and the 1.275 MeV peak from the decay to Ne-22 are visible.

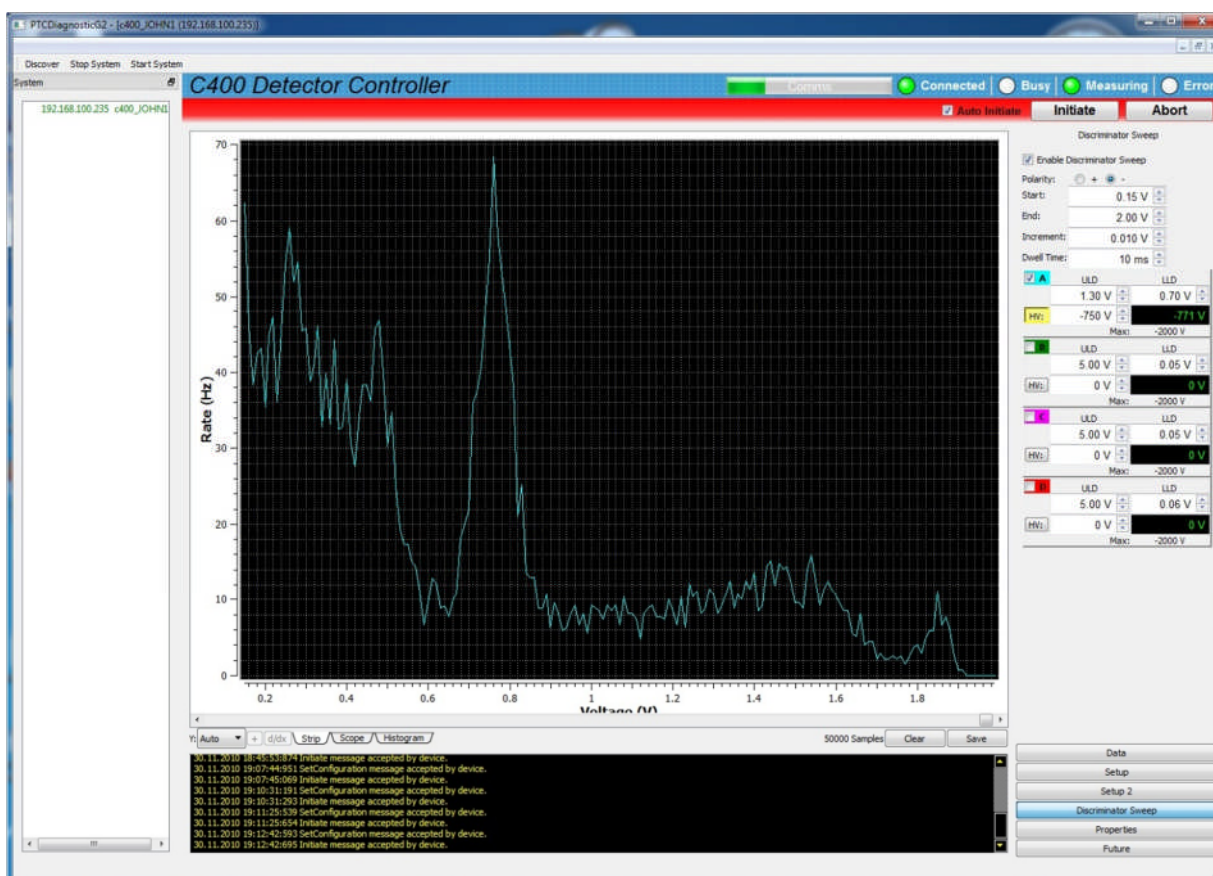


Figure 34. Na-22 spectrum acquire using discriminator sweep.

Four-channel Pulse Counting Detector Controller

There are two unwanted features in pulse height spectra that you should avoid. If you start the sweep at very low voltage, below about 0.05 V, you will start to see electronic noise from the pre-amplifier. It is easy to identify because it will be present even when there is no high voltage bias on the detector. Remove the noise from the spectrum by starting the sweep just above the noise level.

Secondly, be aware of the maximum height that pulses can reach. If you are using a CP10 pre-amplifier, the maximum voltage output is about 2 V. Any larger pulses will saturate at this level, resulting in a spurious peak at the high voltage end of the spectrum. Even if your detector gain is set so that the pulses of interest to you are well clear of saturation, there will still be a natural background of high energy radiation. Avoid creating a spurious peak by ensuring the maximum of the sweep is below the maximum output voltage of the pre-amplifier when driving into the 50 ohm load of the C400 input.

Four-channel Pulse Counting Detector Controller

DEADTIME CORRECTION

MODELS OF DEADTIME

All pulse counting detector systems have an unavoidable deadtime. There are two models of deadtime, called non-paralyzable and paralyzable deadtime. In the non-paralyzable case we assume there is a deadtime when a pulse is being handled. Another pulse arriving during this time is not counted, but has no other impact. In the paralyzable case, we assume that another pulse arriving in the deadtime is not counted, but that it can nevertheless extend the deadtime. The non-paralyzable case is more mathematically tractable, so the extent that we can perform a reliable calculation to estimate the counts that were lost, assuming only that the counting statistics are Poissonian, which is the case for most radiation sources.

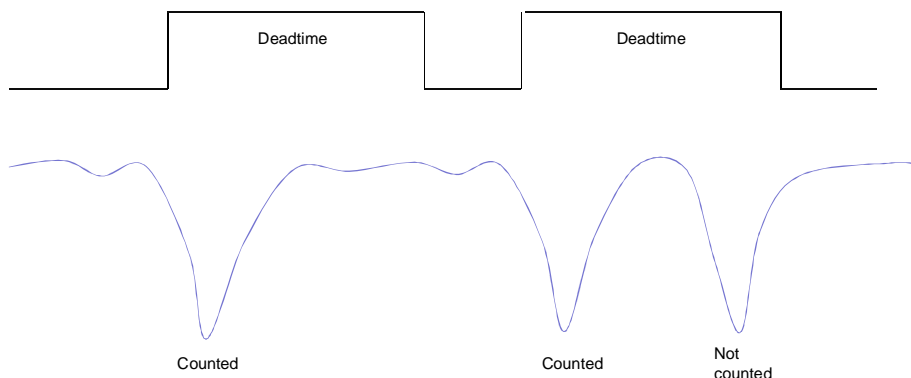


Figure 35. Non-paralysable deadtime model.

The deadtime is composed of individual effects such as the processing time of the electronics and the physical duration of the detector pulse and the lengthening of the detector pulse due to integration on parasitic or deliberately introduced capacitances in the signal processing chain. Generally one component will dominate, and the shorter components have no influence. The deadtime due to the C400 itself is better than 10 nsec, so in most cases the detector/pre-amplifier will be limiting.

The C400 calculates the actual counts for each individual integration period in real time, using the equation

$$N_{\text{actual}} = N_{\text{meas}} / (1 - (\tau/T) * N_{\text{meas}})$$

where N_{meas} is the actual measured counts in integration length T seconds and τ is the fixed deadtime in seconds. The correction can be turned on or off using the check box on the setup page of the PTC Diagnostic program. The correction is generally dependable up to 20 to 30%, but if you are making larger corrections, then it is likely that the non-paralyzable model is becoming invalid as effects like detector gain suppression start to dominate.

Four-channel Pulse Counting Detector Controller

You should exercise care in using data which has been extrapolated too far. As an illustration, if we have a system deadtime of 50 nsec, then an incoming rate of 3.3 MHz will be corrected up to 4 MHz, a 20% correction. If the deadtime is 10 nsec, then the same relative correction occurs at an incoming count rate of 16.7 MHz.

DETERMINING THE DEADTIME

You can often make a reasonable estimate of the deadtime by observing the pulses being delivered to the C400 inputs on a fast oscilloscope, and making an estimate of the typical pulse width. However the most reliable way of measuring the deadtime is to measure the detector response to a radiation source whose intensity you can control at will. Examples include emission stabilized X-ray tubes, or various numbers of small, equal radioactive sources. Alternatively you may be able to control the geometry of the measurement setup to adjust the flux reaching the detector, for example using collimators of known apertures, or by altering the separation between the source and the detector.

The following example uses data from a LaBr_3 scintillation detector measuring the flux from a small stabilized X-ray tube. The geometry and the high voltage on the tube were kept constant, and the tube emission current was adjusted.

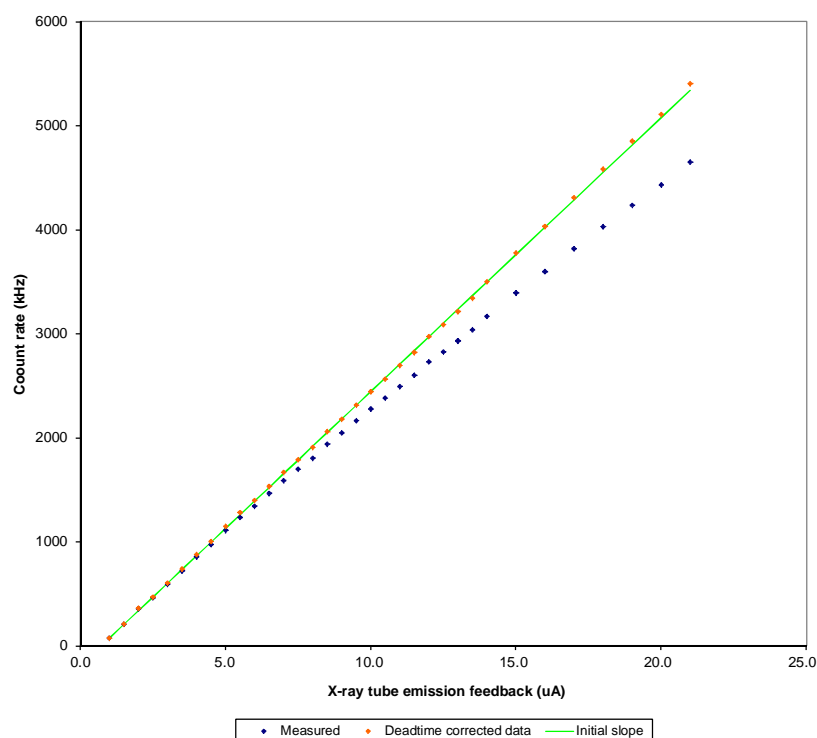


Figure 36. Deadtime measurement using a controllable intensity source.

Four-channel Pulse Counting Detector Controller

The blue points are the measured counts recorded by the C400 with no deadtime correction applied, as a function of the X-ray tube emission current. The green line is a linear fit to the first few data points, where we can assume the effect of deadtime is negligible, and thus represents the ideal detector response. You can see how the recorded counts fall increasingly below the ideal line at increasing intensity, as the losses due to deadtime become significant. The red points are the recorded values with an off-line correction using a deadtime of 30 nsec. This was in accordance with the pulse widths as observed on an oscilloscope. The validity of the proposed value was checked by entering the value as the deadtime in the C400, and repeating the measurements with correction enabled.

Four-channel Pulse Counting Detector Controller

PRE-AMPLIFIER SELECTION

The C400 is compatible with many pre-amplifiers, and can provide power to any that require ± 12 VDC at no more than 100 mA. The Pyramid CP10-A and CP10-B types are designed for use with the C400.

It is important to choose an appropriate type for your application. The C400 is a fast pulse counting device, so an integrating pre-amplifier is generally not suitable, because pulses are converted into steps. A wide-band amplifier will preserve the fast pulse from your detector. This means a low input impedance is required, usually 50 ohms, to minimize unwanted integration of the pulses on the capacitance of the detector plus the cable connecting it to the pre-amplifier. DC-coupled pre-amplifiers will suffer less from baseline shift at high count rates than AC-coupled ones.

You might be inclined to always use the highest bandwidth pre-amplifier available, but this is not the best strategy. The high-speed usually comes with a price in noise levels and DC offset. Most important, however, is that excess bandwidth may simply expose characteristics of the detector output pulses that are detrimental to best pulse counting performance. The following single pulses captured on a digital oscilloscope show the pulses from a LaCl_3 scintillation detector amplified by the CP10-A and the faster CP10-B pre-amplifiers. The bandwidth of the CP10-A is suited to this type of detector in the pulse counting application, but the CP10-B reveals individual electron avalanches in the photomultiplier during the light decay which could confuse the discriminator logic and lead to incorrect counting.

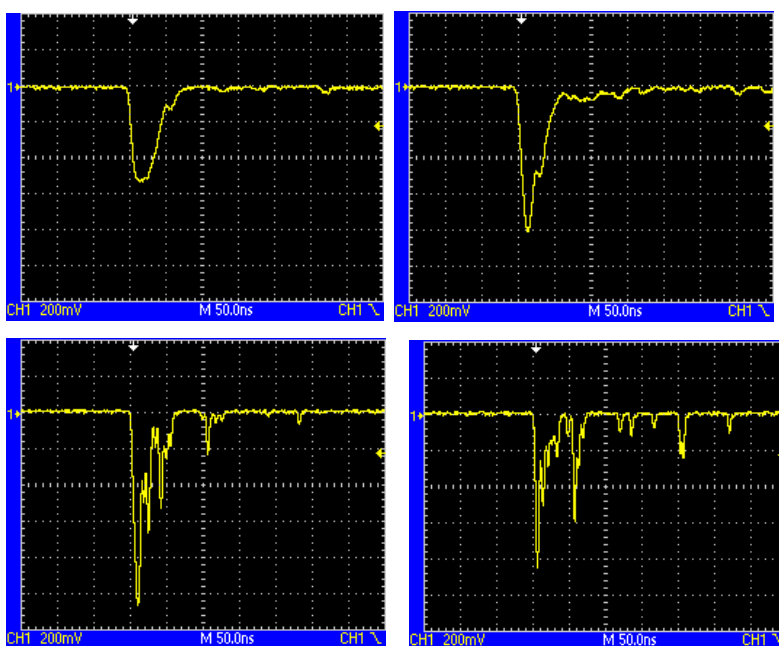


Figure 37. Pulse shapes from a LaCl_3 scintillation detector as seen by the C400 with CP10-A and CP10-B pre-amplifiers.

Four-channel Pulse Counting Detector Controller

PULSER OUTPUTS

The C400 pulser output allows you to perform a full validation of your counting system in suitably configured detector systems. An LED is fitted inside the housing of a scintillation detector so that it can deliver light pulses to the scintillator. If the LED is driven by short pulses, then the resulting detector response can be a good simulation of ionizing radiation events, except with exactly known frequency.

The C400 can generate 5V pulses as short as 1 nsec with a frequency of up to 1 MHz, and current up to 30 mA. The size of the detector response will depend on the pulse duration and the LED current. You can set the current by choosing a series resistor for the LED. For a typical LED forward voltage of 2V, the LED current will be $3000/R$ mA. Thus a series resistor R of 150 ohms will give current of 20 mA, and this is also a good termination for a twisted pair line. However you can try any value within the current compliance of the C400 and the rating of the LED in the detector.

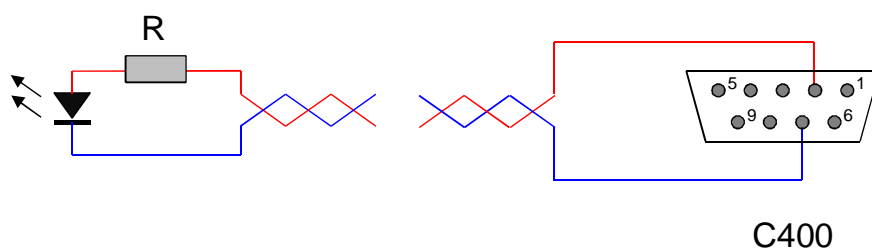


Figure 38. Test pulser LED connection.

Having chosen a suitable LED current, you should then observe the detector pulses on an oscilloscope and adjust the pulse duration to give the best match to real radiation events.

Four-channel Pulse Counting Detector Controller

ENCODER INPUT

Many situations require that the measurements from one or more detectors are correlated with the physical position of a sample, or of the detector itself. The encoder input on the C400 allows a standard position encoder to be read out in synchronization with the pulse counting measurements.

The encoder input handles standard differential quadrature signals.

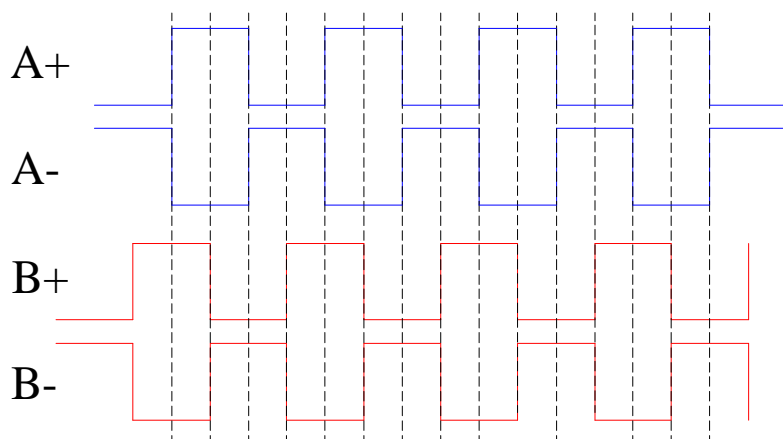


Figure 39. Differential quadrature signals from a shaft encoder, clockwise movement.

The input voltages can be up to ± 25 V, with a differential of greater than 0.2 V necessary to define an input state. Input levels of 0 V and 5 V are typical. A +5 VDC output is provided which can be used to power suitably-rated encoders. A differential index signal is also supported.

Differential signals are recommended for their good noise immunity. However many encoders only provide A and B inputs. In this case you should bias the unused A- and B- inputs (and Z- if you are using the index) to a voltage halfway between the high and low states of the pulses.

Four-channel Pulse Counting Detector Controller

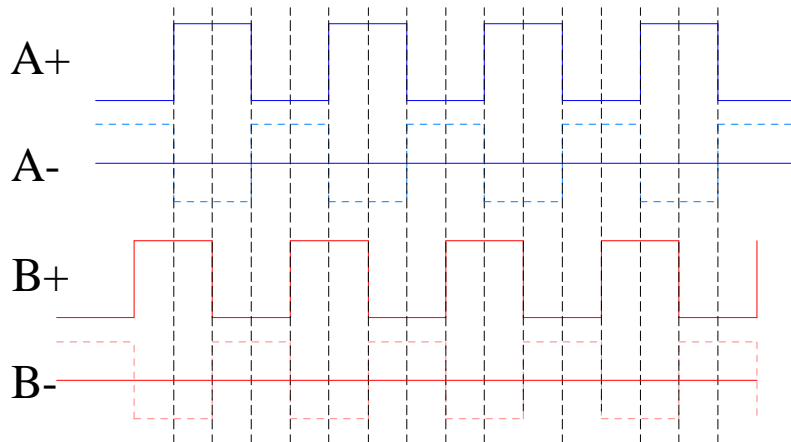


Figure 40. Encoder signals when no differential signals are available.

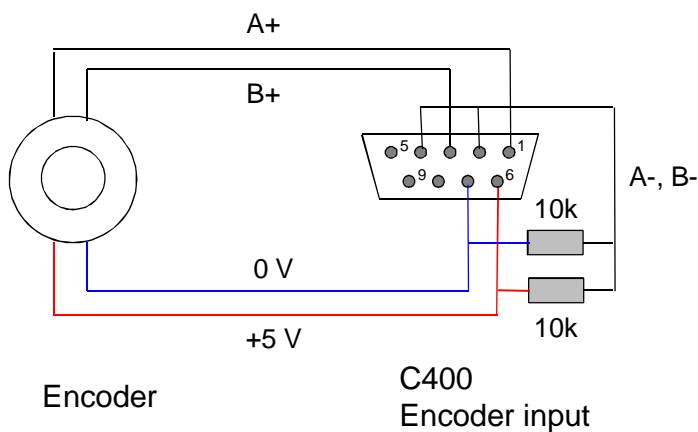


Figure 41. Example bias circuit for encoder inputs when no differential signals are available.

Four-channel Pulse Counting Detector Controller


CONNECTORS

FRONT PANEL CONNECTORS

Four identical sets of connectors for the four channels.


Signal input

Four Lemo 00 50 ohm coaxial.

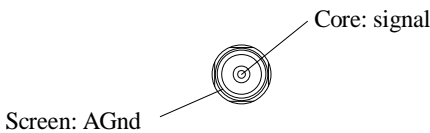

Pre-amplifier power and other services


Four Dsub 9 pin female.

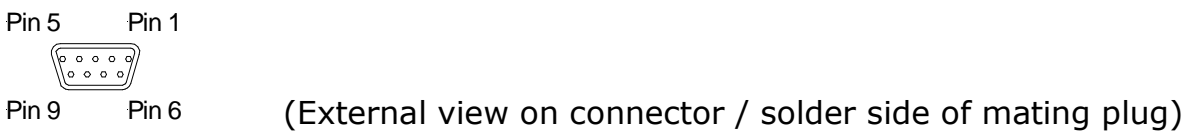

1	DGnd	6	DGnd
2	Pulser out	7	DGnd
3	AGnd	8	+12 VDC
4	-12 VDC	9	AGnd
5	Analog output		

High voltage outputs

Four SHV male. To mate with standard SHV connector.

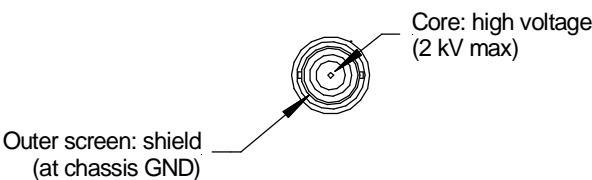

REAR PANEL CONNECTORS

Ethernet communications

RJ-45 jack. To mate with standard RJ-45 plug.

Auto MDIX facility - cable can be direct or crossover type.

Four-channel Pulse Counting Detector Controller


Fiber-optic communications:

HFBT ST bayonets suitable for 1 mm plastic or 200 μ m silica fiber. 664 nm (visible red) light.

Dark casing = receiver, light casing = transmitter.

RS-232 / RS-485 communications

Six pin mini-DIN socket (PS/2 mouse/keyboard type).

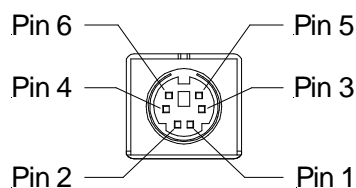

(External view on connector / solder side of mating plug)


1	RS-232 Tx / RS-485 Tx-	4	n/c
2	RS-232 Rx / RS-485 Rx+	5	RS-485 Tx+
3	Gnd	6	RS-485 Rx-

The socket incorporates a sensor switch that allows the C400 to detect that a plug has been connected. When a connection is made, the RS-232 / RS-485 transceiver is active, and the communication mode is set by the mode switch.

(External view on connector / solder side of mating plug)


Gate output

BNC socket (female). To mate with standard signal BNC. Able to drive a 50 ohm load.

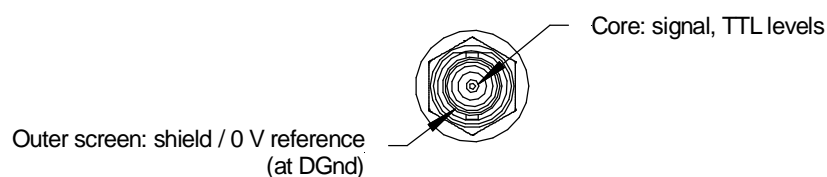

Gate input

BNC socket (female). To mate with standard signal BNC.

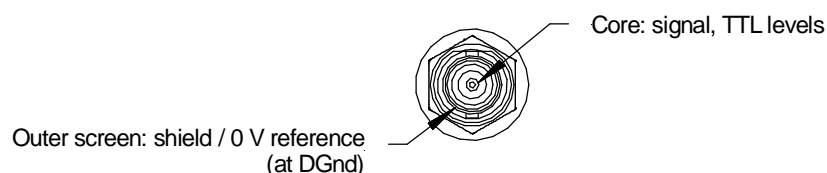

Four-channel Pulse Counting Detector Controller

Encoder input

Dsub 9 pin female.

Pin 5 Pin 1




Pin 9 Pin 6

(External view on connector / solder side of mating plug)

1	Enc A+	6	+5 VDC out
2	Enc A-	7	DGnd
3	Enc B+	8	DGnd
4	Enc B-	9	Enc Z+
5	Enc Z-		


Monitor outputs

Four Lemo 00 coaxial.

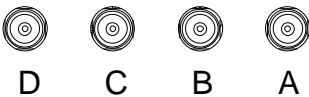

Power input

2.1 mm threaded jack. To mate with Switchcraft S761K or equivalent

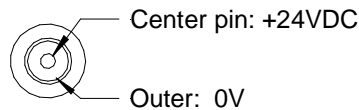

Ground lug

M4 threaded stud. To mate with M4 ring lug.

Four-channel Pulse Counting Detector Controller

CONTROLS AND INDICATORS

FRONT PANEL CONTROLS

None.

REAR PANEL CONTROLS

Reset button

Momentary push-button that forces a warm reset of the on-board processor.

Address switch

16 position rotary switch setting device address for use when you have the C400 on a fiber optic loop. Choice of address is arbitrary, but each device in a fiber-optic loop system must have a unique address.

<i>Setting</i>	<i>Function</i>
0	(Reserved to loop controller)
1-15	Available address settings.

Note: Fiber optic loop communication will be added in a future software release.

Mode switch

10 position rotary switch setting communications mode. This switch is read if a connection is made to the serial connector. The ASCII protocol is provided for ease of connection to existing systems and simple terminal programs.

<i>Setting</i>	<i>Function</i>
0	8 bit binary, 115 kbps, RS-232
1	8 bit binary, 57.6 kbps, RS-232
2	8 bit binary, 19.2 kbps, RS-232
3	ASCII, 115.2 kbps , RS-232
4	ASCII, 57.6 kbps , RS-232
5	ASCII, 19.2 kbps, RS-232
6	8 bit binary, 115 kbps, RS-485
7	8 bit binary, 57.6 kbps, RS-485
8	ASCII, 115.2 kbps , RS-485
9	ASCII, 57.6 kbps , RS-485

Four-channel Pulse Counting Detector Controller

FRONT PANEL INDICATORS

HV on

Four yellow LEDs. Illuminated if the HV supply is enabled on the channel.

REAR PANEL INDICATORS

Eight green LEDs.

24V   Com

2.5V   Ethernet

Initiated   Optical

Active   Serial

24V

Green LED. 24 VDC power is present.

2.5V

Green LED. 5V DC-DC converter is running, 2.5 V regulator running and generating 2.5 VDC power.

Initiated

Green LED. Unit has been initiated and triggered to acquire data.

Active

Green LED. C400 processor has booted and is running normally.

Com

Green LED. A host communication channel is active.

Ethernet

Green LED. Ethernet commands are being received. The on state is stretched so that it can be seen. This LED does not illuminate for pings or if the C400 is streaming data to the host.

Optical

Green LED. Fiber optic is the active communication interface.

Serial

Green LED. Serial (RS-232 or RS-485) is the active communication interface.

Four-channel Pulse Counting Detector Controller

Power up sequence

On power-up, all LEDs light initially. The Initiated and Active LEDs then extinguish, and the Com, Ethernet, Fiber and Serial LEDs light in sequence while the C400 is booting.

RJ-45 indicators

Green LEDs incorporated in the Ethernet connector show that a connection has been established (left hand LED) and that messages are passing (right hand LED).

Four-channel Pulse Counting Detector Controller

COMMUNICATIONS INTERFACES

The C400 is a member of the PSI G2 range of devices. The unit is provided with four hardware interfaces, RS-232, RS-485, fiber-optic and Ethernet. The RS-232 and RS-485 interfaces are intended for simple direct connection to PCs, with no other equipment necessary. The fiber-optic interface provides greater speed, excellent noise immunity, and allows multiple devices to be connected in a looped topology. It requires a fiber-optic adaptor or loop controller device to connect to the host computer. The fiber-optic interface is well-suited to large systems and experiments. The Ethernet interface is provided to allow direct connection to systems which are integrated using standard local area network hardware and protocols.

The C400 can support multiple clients, and in particular you can be connected via Ethernet and the serial port at the same time.

If the serial connector is attached, then the serial port becomes active, irrespective of any other connections. The baud rate, type of protocol and choice of RS-232 or RS-485 levels is made using the mode switch.

If the Ethernet connection is made then the Ethernet port becomes active. The mode switch and address switch have no function when using Ethernet communications. LEDs incorporated in the RJ-45 jack indicate activity on the port.

If neither serial nor Ethernet connections are made, then the fiber optic port is the active interface. The mode switch has no function. The address switch sets the fiber optic loop address.

RS-232 CONNECTION

A basic three-wire RS-232 connection is all that is required to communicate with the C400. The C400 uses a 6-pin mini-DIN connector. An adaptor to a 9-pin female DSub is available from Pyramid; a standard pin to pin male to female 9-pin DSub cable can then be used to connect to the PC. The cable length should not exceed about 5m, especially at the maximum baud rate. For longer serial cable runs you should use RS-485.

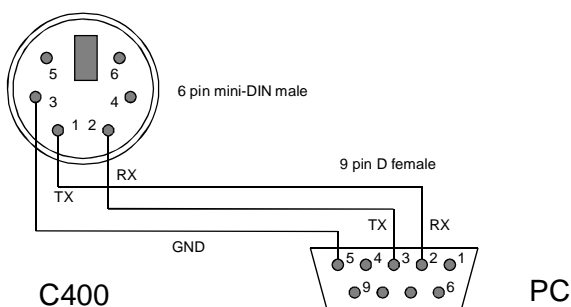


Figure 42. RS-232 cable C400 to PC.

Four-channel Pulse Counting Detector Controller

RS-485 CONNECTION EXAMPLE

RS-485 is used to extend a serial connection over a long distance. If you have a RS-232 serial port on the host PC, then an RS-232 to RS-485 converter is used at the host PC end. The built-in RS-485 of the C400 allows a direct connection at its end.

A commonly-used converter is the MOXA TC100. The Moxa TCC-80 is a low-cost alternative for less critical applications. The converter should be configured for four wire (full duplex) RS-485 operation. The C400 provides parallel termination for transmit and receive. It is therefore optional whether you also terminate at the TC100 end, but there is no harm if you do. The recommended Dip switch setting for the TC100 is:

Sw1	Sw2	Sw3
OFF	ON	ON

The recommended Dip switch setting for the TCC-80 is:

Sw1	Sw2	Sw3
ON	OFF	ON

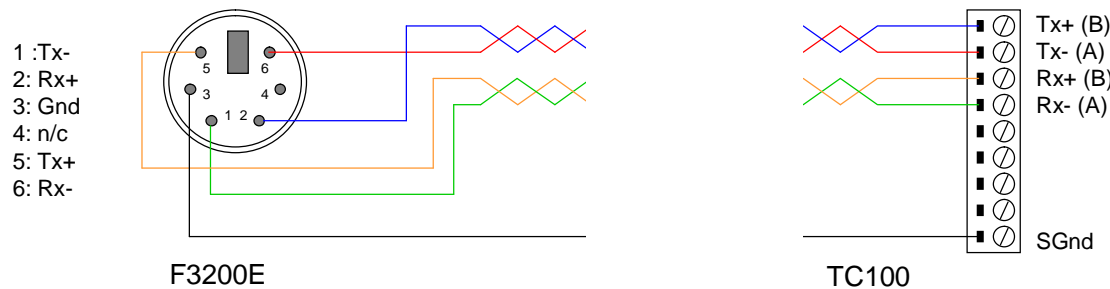


Figure 43. RS-485 cable C400 to TC100.

In order to enable the RS-485 interface on the C400, you must have a physical connection to the serial port, and you must select one of the relevant modes (modes 6 through 9). If you have trouble getting the RS-485 connection to work, check carefully to ensure that you have all the transmit and receive lines connected as shown.

Four-channel Pulse Counting Detector Controller

ASCII COMMUNICATION

The C400 is a high performance device that will give you best performance with a dedicated host software application and high bandwidth communications interface. Nevertheless, there are numerous applications that do not require especially high performance, and connections to the C400 may be most convenient with a simple serial interface. The C400 supports a basic ASCII interface over RS-232 or RS-485. It supports some of the commands prescribed by Standard Commands for Programmable Instruments (SCPI) and IEEE 488.2, plus specific commands as required by the operation of the device.

All commands to the C400 should be terminated with LF only. All C400 responses are terminated with CR LF. The backspace is CTRL-H (ASCII char 8).

The commands are grouped with a hierarchical structure, with the levels separated by the colon character. For example:

CONFigure:PERiod 1e-2// This command configures the integration period to have a length of 10 milliseconds.

SCPI provides for a long and short form for each command. The short forms are indicated by the capitalized part of the command.

Commands which have a query equivalent for readback are marked with "(?)" in the following tables. Parameters are generally passed to the C400 with the set version of the command, but no parameters are passed for the query version. For example:

CONFigure:PERiod 1e-2// set the integration period to 10 msec
CONFigure:PERiod? // query the integration period

In the following tables,

{ } denotes a required argument.

[] denotes an optional argument.

Four-channel Pulse Counting Detector Controller

IEEE 488.2 MANDATORY COMMANDS

*CLS		Not supported
*ESE	(?)	Not supported
*ESR		Not supported
*IDN?		Identification Query. C400 returns manufacturer, model number, serial number, firmware version
*OPC	(?)	Not supported
*RST		Reset Command. Restart the device software and return to the *RST default conditions. This command is not yet supported.
*SRE	(?)	Not supported.
*STB?		Not supported
*TST?		Not supported
*WAI		Not supported

IEEE 488.2 OPTIONAL COMMANDS

*RCL		Recall device configuration from flash
*SAV		Save device configuration to flash

The settings covered by *RCL and *SAV are:

CONFigure:DAC

CONFigure:DLO

CONFigure:DHI

CONFigure:POLarity

CONFigure:PERiod

CONFigure:PULSer

CONFigure:DEADtime

CONFigure:HIVoltage:VOLts (Note that the CONFigure:HIVoltage:ENable settings are not included).

TRIGger:BUFFer

TRIGger:BURst

TRIGger:MODE

TRIGger:SOURce

TRIGger:POLarity

Settings saved using *SAV are automatically recalled when the device is turned on.

Four-channel Pulse Counting Detector Controller

C400 SPECIFIC COMMANDS

ABORt				Abort measurement
CONFigur e	:DAC {Ch} {Val}		(?)	Set analog outputs DAC Ch to Val volts. (Query all 4 DAC settings.).
	:DEADtime		(?)	Set (query) the deadtime correction in ns to be applied to the incoming counts for all channels.
	:DHI {Hi1 Hi2 Hi3 Hi4}		(?)	Set (query) the discriminator high level Hi1..Hi4 in volts for each channel. The sign is ignored (see CONFigure:POLarity).
	:DLO {Lo1 Lo2 Lo3 Lo4}		(?)	Set (query) the discriminator low level Lo1..Lo4 in volts for each channel. The sign is ignored (see CONFigure:POLarity).
	:ENCODer {Enc}			(Command not yet supported)
	:HIVoltage	:MAXvalue {HVm1 HVm2 HVm3 HVm4}	(?)	Set (query) maximum allowable external high voltage setting HVm1..HVm4 in volts for each channel. This is a software enforced limit.
		:SUPply?		Query maximum allowable external high voltage setting in volts for each channel, as determined by the installed HV module.
		:VOLts {HV1 HV2 HV3 HV4}	(?)	Set (query) the auxiliary (external) high voltage setting HV1..HV4 in volts for each channel. See FETch:HIVoltage? for readback values.
		:ENable {0 1} {0 1} {0 1} {0 1}	(?)	Enable/disable (query) the auxiliary (external) high voltage power supply for each channel.
	:PERiod {Per}		(?)	Set (query) integration period Per in seconds.
	:POLarity {N P} {N P} {N P} {N P}		(?)	Set (query) the discriminator polarity for each channel.
	:PULser {Period Width}		(?)	Set (query) the LED pulser output. Period and width are in nsec.
FETch	:COUNTs? [numSample s]			Returns counts for the numSamples most recent triggers for the four channels. Returns: numSamples *{integration time,count1,count2,count3,count4,timesta mp,trigger count,Lo1,Lo2,Lo3,Lo4 CR LF}. It is only possible to request multiple samples if TRIGger:BUFFer > 0. Otherwise, numSamples should be omitted.

Four-channel Pulse Counting Detector Controller

	:RATE? [numSamples]			Returns rates for the numSamples most recent triggers for the four channels. Returns: numSamples * {integration time, rate1, rate2, rate3, rate4, timestamp, trigger count, Lo1, Lo2, Lo3, Lo4 CR LF}. It is only possible to request multiple samples if TRIGGER:BUFFER > 0. Otherwise, numSamples should be omitted.
	:DIGITAL?			Read digital status bits. See following section for the meaning of the bits.
	:ENCODER?			Return the current encoder count.
	:HIVOLTAGE?			Read HV readback ADCs. Returns {HV1, HV2, HV3, HV4}.
INITIATE				Initiate data acquisition on valid trigger.
PULSE {channels}				Send out LED pulses on channels (0-15 binary pattern).
SCAN {start_V} {stop_V} {window_V} {period}				Configure and initiate a discriminator sweep starting at start_V (signed) and ending at stop_V (signed) in increments of window_V and a dwell time of period seconds.
SYSEM	:COMMUNICATION	:TIMEOUT	(?)	Set (query) timeout in seconds. 0 = timeout disabled. C400 will turn HV off if no valid message is received in the timeout period. This command is not yet supported.
	:SERIALNUMBER?			Query the serial number of the C400.
	:VERSION?			Query the SCPI standard version.
TRIGGER	:BUFFER {size}		(?)	Set (query) the C400 on-boards data buffer size.
	:BURST {count}		(?)	Set (query) the number of samples acquired in a burst (i.e. before a pause). This is 1 in EXTERNAL_START_HOLD mode.
	:MODE {mode}		(?)	Set (query) the trigger mode to mode. The options are: CUSTOM, INTERNAL, EXTERNAL_START, EXTERNAL_START_STOP, EXTERNAL_START_HOLD, EXTERNAL_WINDOWED, DISCRIMINATOR_SWEEP
	:POLARITY {0 1}		(?)	Set (query) external gate polarity (external trigger only). 0 = rising edge, 1 = falling edge (invert BNC trigger)
	:SOURCE	:START {INTERNAL B	(?)	Set (query) the start trigger source. Used in CUSTOM, EXTERNAL_START,

Four-channel Pulse Counting Detector Controller

		NC})	EXTERNAL_STARTSTOP, EXTERNAL_START_HOLD, and EXTERNAL_WINDOWED modes.
		:STOP {INTERNAL B NC}	(?)	Set (query) the stop trigger source. Used in CUSTOM and EXTERNAL_STARTSTOP modes.
			(?)	Set (query) the pause/burst trigger source. Used in CUSTOM and EXTERNAL_WINDOWED modes.

DIGITAL STATUS BITS

0	Connected
1	Pending/Busy
2	Error
3	Reserved
4	Output changed (NOT SUPPORTED)
5	Reserved
6	Urgent device error (NOT SUPPORTED)
7	Reserved
8	Updating firmware
9	Starting up
10-15	Reserved
16	Measuring

USING PuTTY TO TEST THE SERIAL ASCII CONNECTION

Microsoft Windows no longer includes a terminal emulator that you can use to test the ASCII communications. Various public domain programs are available to perform this task. We have tested PuTTY and Realterm with the C400. The following screenshots show the configuration of PuTTY with 115kbps baud rate.

Four-channel Pulse Counting Detector Controller

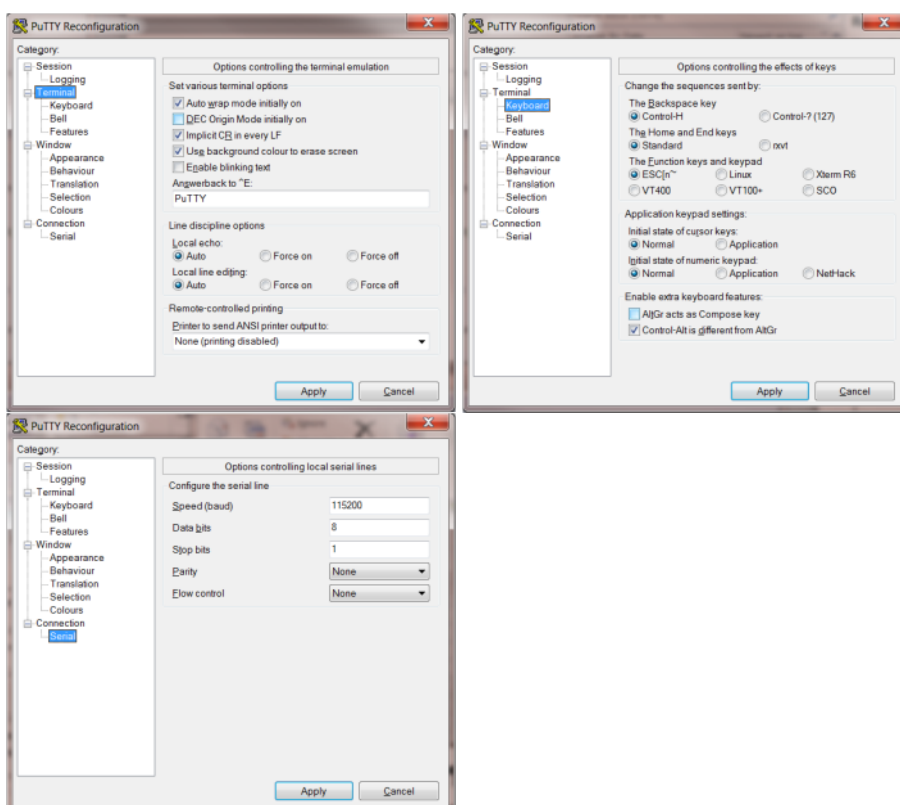


Figure 44. Setup of PuTTY for ASCII serial communication with the C400.

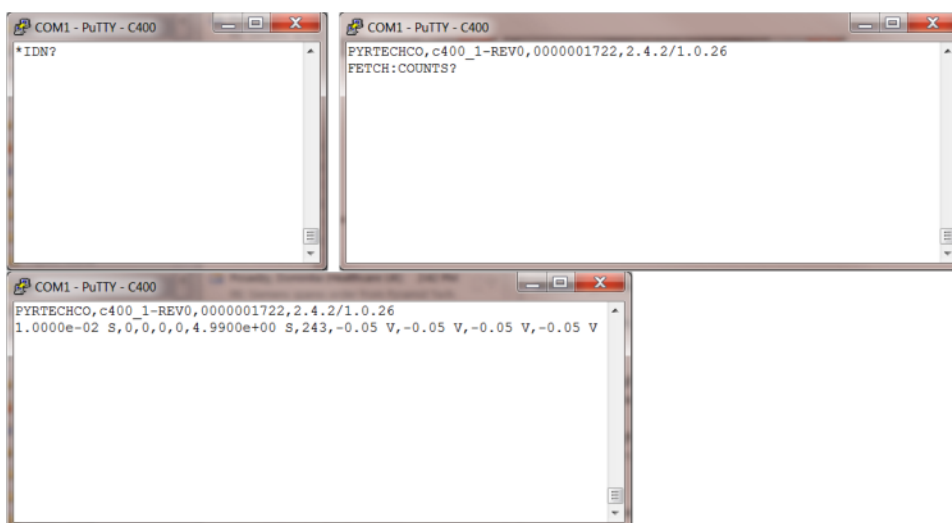


Figure 45. Example terminal session.

Four-channel Pulse Counting Detector Controller

CREATING A HOST SOFTWARE APPLICATION

Pyramid Technical Consultants, Inc. provides a full application programming interface for the C400, based upon a C++ library. Full details are available upon request.

Four-channel Pulse Counting Detector Controller

FAULT-FINDING

Symptom	Possible Cause	Confirmation	Solution
High background count rates	Discriminator lower level too low	Increase discriminator – counts drop suddenly	Set low discriminator limit to exclude electronic noise.
	Noisy pre-amplifier	Change pre-amplifier	Set low discriminator limit to exclude electronic noise for pre-amp in use.
	Large DC offset on pre-amplifier	Change pre-amplifier	Set low discriminator limit to exclude electronic noise for pre-amp in use.
No pulses	Pulse polarity set incorrectly	Look at pulses on oscilloscope (50 ohm terminated)	Use correct polarity setting.
	Pulses are outside discriminator window	Open window settings and recheck. Look at pulses on oscilloscope (50 ohm terminated)	Use appropriate discriminator settings. Discriminator sweep function can assist you to select the values.
Counts higher than expected	Deadtime correction is enabled; deadtime setting so too high.	Disable deadtime correction.	Ensure use of deadtime correction is appropriate and the correct value is entered.
	Pulser output is connected and enabled.	Turn off pulser.	Disable pulser when taking data.

Four-channel Pulse Counting Detector Controller

Counts lower than expected.	Detector gain has dropped.	Look at pulses on oscilloscope (50 ohm terminated). Increase high voltage and/or reduce discriminator lower level.	Use correct high voltage setting. Replace detector if needed.
No signal on monitor outputs.	Unit is not measuring.	Signal appears when unit initiated and triggered.	Ensure C400 is measuring when you want to see the monitor output.
Unit not collecting data	Inappropriate trigger mode	Select internal trigger mode.	Use appropriate trigger mode.
Encoder not reading, or erratic.	Encoder is not producing differential pulses, and A-, B- lines are floating.	Check pulses with oscilloscope.	Use a differential encoder, or tie the unused lines to a voltage halfway between the two input states.
Encoder counting backwards	Direction of motion can be defined arbitrarily for incremental encoders.		Check invert encoder option (if using PTC Diagnostic).
Encoder counts not updating	C400 not initiated	Initiate C400 and check encoder counts.	
	Communication link timeout		Investigate and fix communications issue. Use a longer timeout setting.
No or incorrect response to external gate	Incorrect gate polarity selected.		Use correct polarity.

Four-channel Pulse Counting Detector Controller

No or low high voltage	Shorted to ground in external circuit	Monitor HV reading zero or very low relative to setpoint. Monitor value recovers if C400 disconnected from the external circuit.	Eliminate shorts to ground.
	External load resistance too low.	Voltage recovers if setting is reduced.	Do not attempt to draw more than 1 watt from each high voltage supply. Higher power modules are available to special order.
High voltage not at setpoint	A high compliance source such as a charged particle beam is driving the HV electrode.	Monitor value recovers if C400 disconnected from the external circuit.	Change geometry to reduce beam strike.
Cannot set high voltage	Trying to set above the maximum allowed value soft limit.	Sets OK if a lower value is chosen.	If allowed, increase the maximum allowed value.
Unable to communicate via Ethernet	Incorrect IP address for C400 or host (not in the same subdomain).	Check settings of C400 host PC.	Use consistent IP addresses.
	Messages being blocked by anti virus software.	Disable anti-virus software	Set up allowed channels for C400 messages.
Unexpected changes to C400 state	Another host is communicating with the C400.	Change IP address. Use a direct cable connection instead of a network.	Set up IP addresses and subnet masks to prevent conflicts.

Four-channel Pulse Counting Detector Controller

Communications interruptions	Other processes on PC host interfering with comms ports.		Use a dedicated PC with simple configuration and minimum number of processes running.
Unable to connect on serial port	Another program is using the COM port.	Try to access the required port with a terminal program.	Choose another port or close down the other program.
	Incorrect port settings.		Correct the settings.
	Incorrect cable.		Make up a suitable cable.
Unexpected behaviour	Unit is overheating due to failed fan	Unit feels warm. Fan stopped.	Return unit for fan replacement.
	Unit is overheating due to airflow blockage	Unit feels warm. Performance improves when operated on the bench.	Clear space around the unit and improve airflow.

Four-channel Pulse Counting Detector Controller

MAINTENANCE

The C400 does not require routine maintenance, except to clear any dust accumulation in the fan filter.

There are no user-serviceable parts inside.



CAUTION. High voltages may be present inside the case. Do not open the case when power is applied.

If there is buildup of dust in the filter, you should clear this by vacuum cleaning in situ, or by removing the filter element and cleaning it separately with an air jet. Note that detaching the filter element also detaches the fan from the case.

The C400 is fitted with a 1.1 amp automatically resetting positive temperature coefficient (PTC) fuse in the 24 VDC input. There are individual 200 mA PTC fuses on the pre-amplifier power outputs and the fan power feed. No user intervention is required if the fuse operates due to overcurrent. The fuse will reset when the overcurrent condition ends.

Four-channel Pulse Counting Detector Controller

RETURNS PROCEDURE

Damaged or faulty units cannot be returned unless a Returns Material Authorization (RMA) number has been issued by FMB Oxford. If you need to return a unit, contact enquiries@fmb-oxford.com, stating:

- model
- serial number
- nature of fault

An RMA will be issued, including details of which service centre to return the unit to.

Four-channel Pulse Counting Detector Controller

SUPPORT

Manual and software driver updates are available for download from the Pyramid Technical Consultants website at www.ptcusa.com. Technical support is available by email from enquiries@fmb-oxford.com or support@ptcusa.com. Please provide the model number and serial number of your unit, plus relevant details of your application.

Four-channel Pulse Counting Detector Controller

DISPOSAL

We hope that the C400 gives you long and reliable service. The C400 is manufactured to be compliance with the European Union RoHS Directive 2002/95/EC, and as such should not present any health hazard. Nevertheless, when your C400 has reached the end of its working life, you must dispose of it in accordance with local regulations in force. If you are disposing of the product in the European Union, this includes compliance with the Waste Electrical and Electronic Equipment Directive (WEEE) 2002/96/EC. Please contact FMB Oxford for instructions when you wish to dispose of the device.

Four-channel Pulse Counting Detector Controller

DECLARATION OF CONFORMITY

(certificate in preparation)